


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Ifosfamide

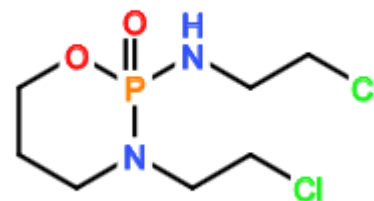

CHROMATOGRAM

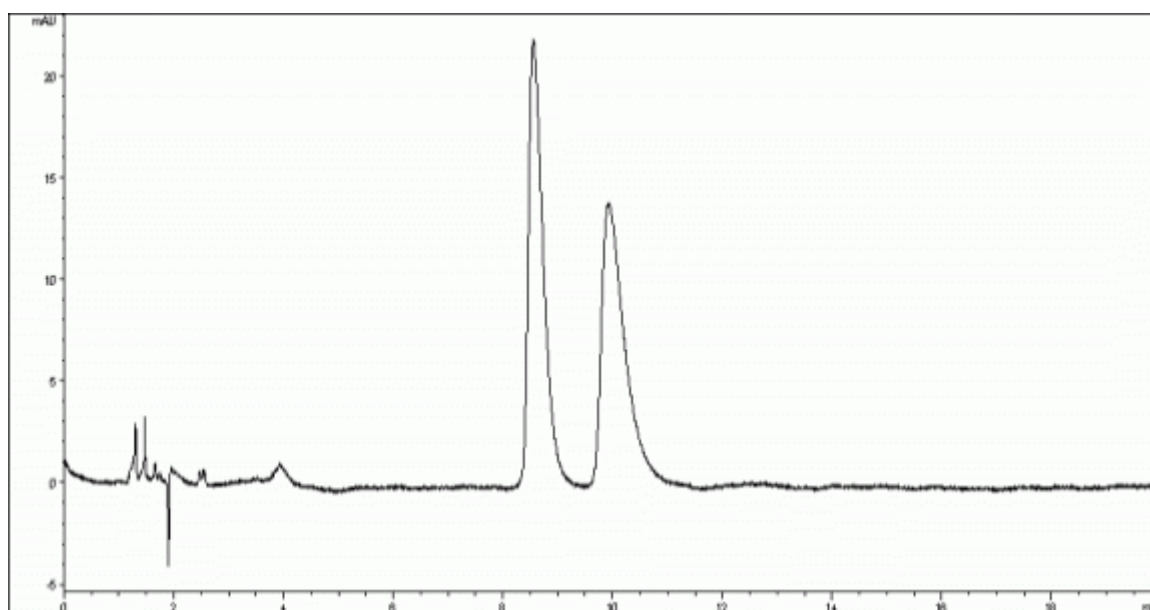

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IA-3		
Column Size (Particle Size)	4.6 x 150 mm (3 µm)		
Mobile Phase Composition	water / acetonitrile = 80 / 20		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 200 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	8.6	3.90	1 - 2	2.21
2	9.9	4.67		