


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Licarbazepine

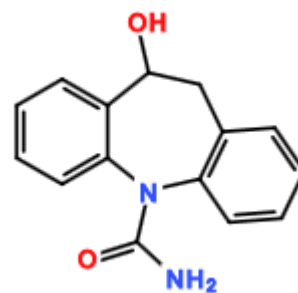

## CHROMATOGRAM

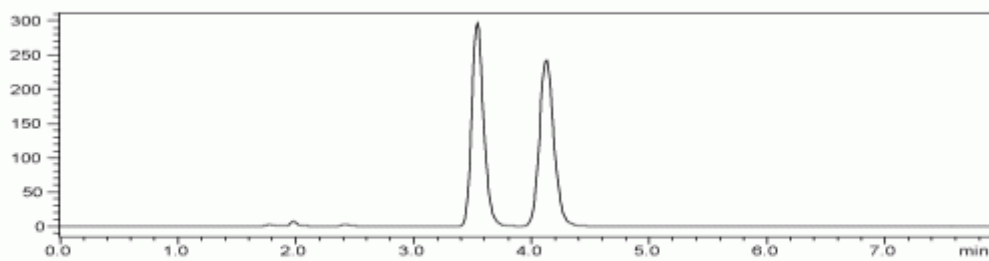

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IC-3		
Column Size (Particle Size)	4.6 x 150 mm (3 µm)		
Mobile Phase Composition	5mM NH4OAc aq. / methanol = 10 / 90		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 220 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	3.5	0.96	1 - 2	2.63
2	4.1	1.28		