


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Omeprazole

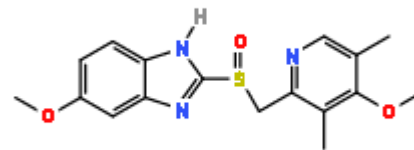

CHROMATOGRAM

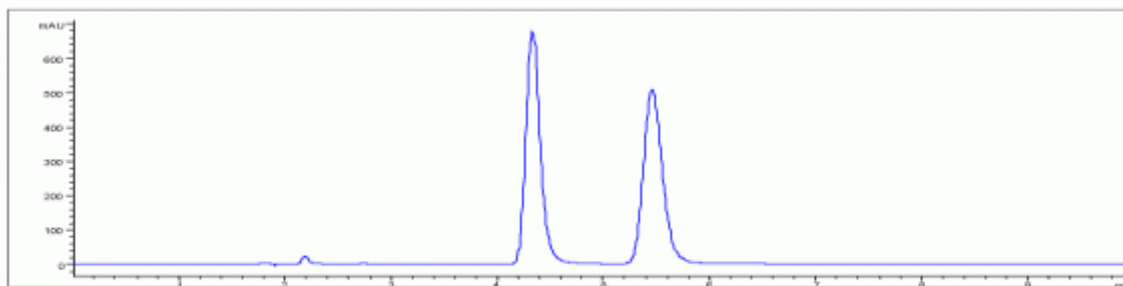

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IC-3		
Column Size (Particle Size)	4.6 x 150 mm (3 µm)		
Mobile Phase Composition	5mM NH ₄ HCO ₃ aq. / methanol = 10 / 90		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 300 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	4.3	1.40	1 - 2	3.53
2	5.5	2.06		