


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

N-CBZ-DL-Serine

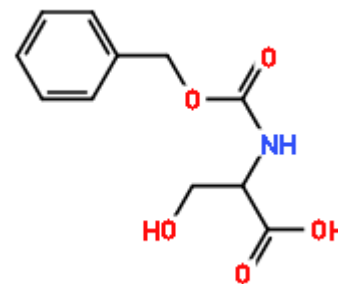

CHROMATOGRAM

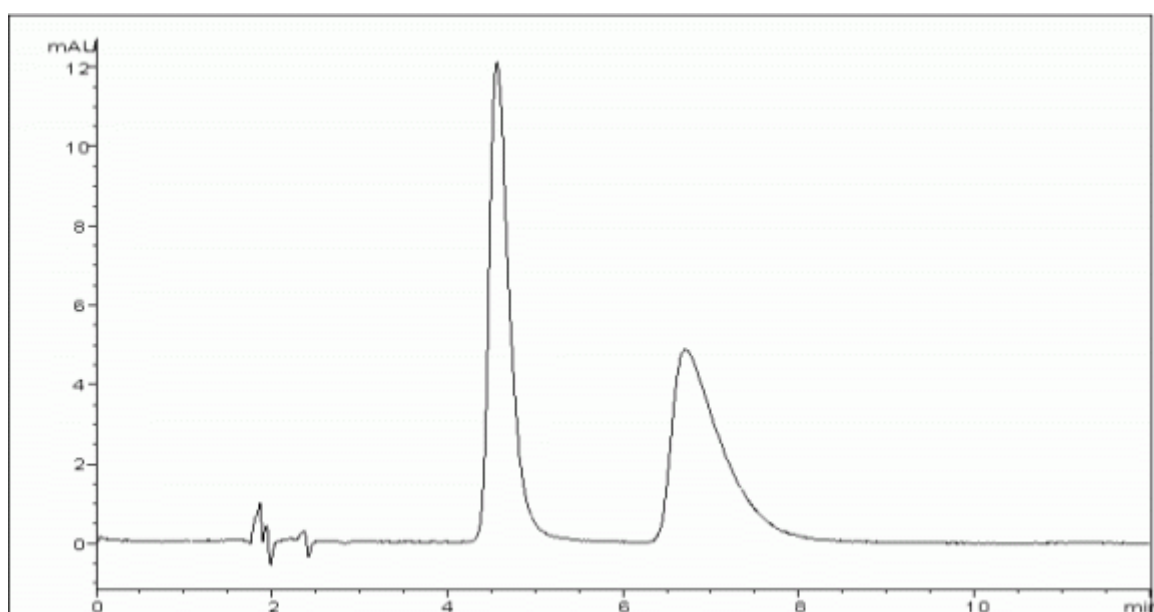

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IA		
Column Size (Particle Size)	4.6 x 150 mm (5 µm)		
Mobile Phase Composition	formic acid aq. pH2.0 / methanol = 55 / 45		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 270 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	4.6	1.31	1 - 2	3.07
2	6.7	2.40		