


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

N-FMOC-DL-Phenylalanine

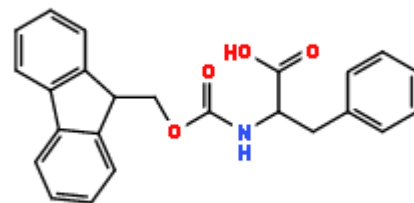

CHROMATOGRAM

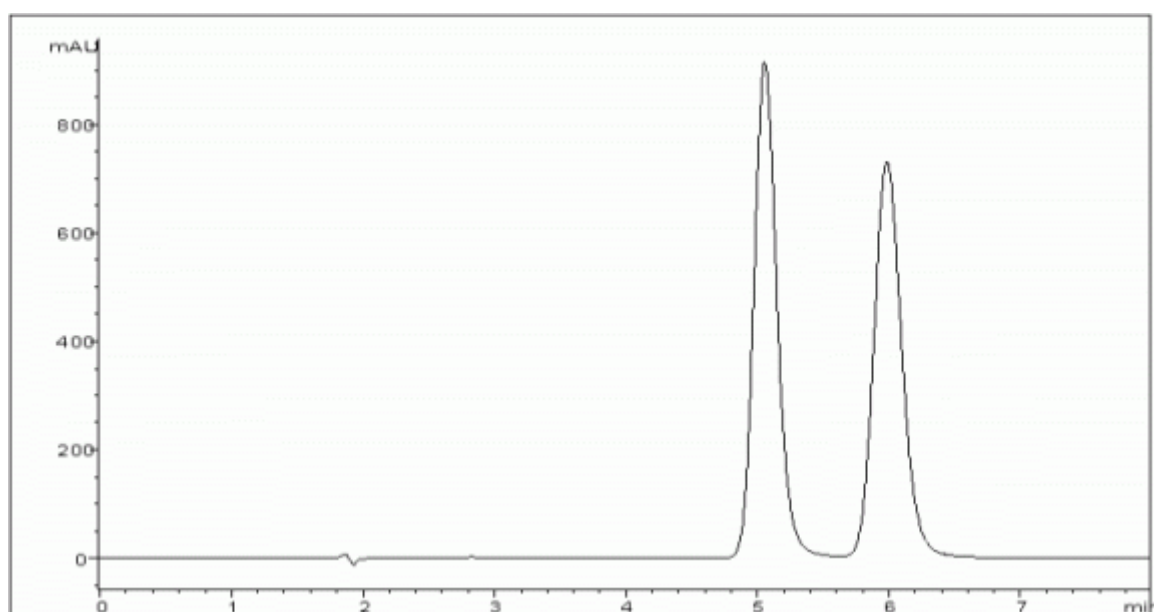

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IC		
Column Size (Particle Size)	4.6 x 150 mm (5 µm)		
Mobile Phase Composition	formic acid aq. pH2.0 / acetonitrile = 50 / 50		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	5.1	1.54	1 - 2	2.68
2	6.0	2.00		