


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Praziquantel

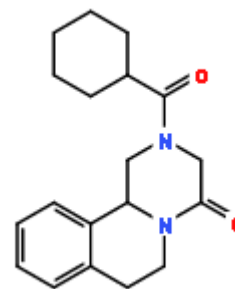

CHROMATOGRAM

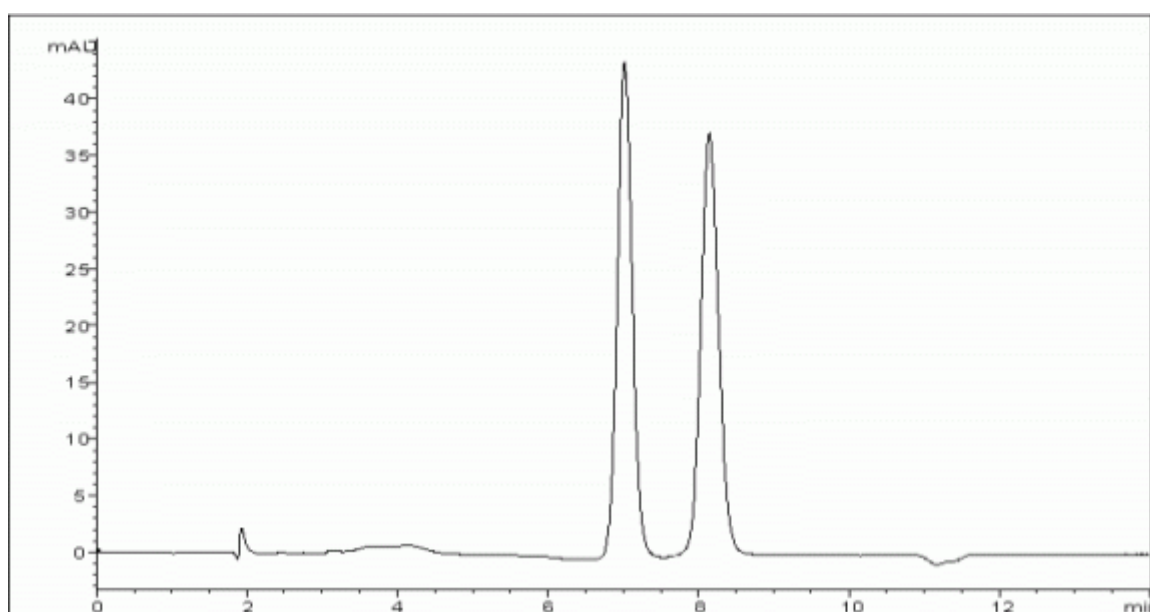

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IB		
Column Size (Particle Size)	4.6 x 150 mm (5 µm)		
Mobile Phase Composition	water / tetrahydrofuran = 65 / 35		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.02 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	7.0	2.60	1 - 2	2.87
2	8.2	3.17		