


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

5-Methyl-5-phenylhydantoin

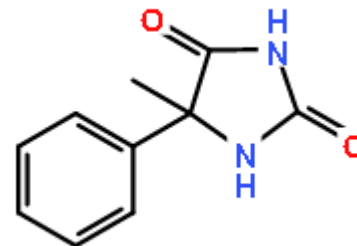

CHROMATOGRAM

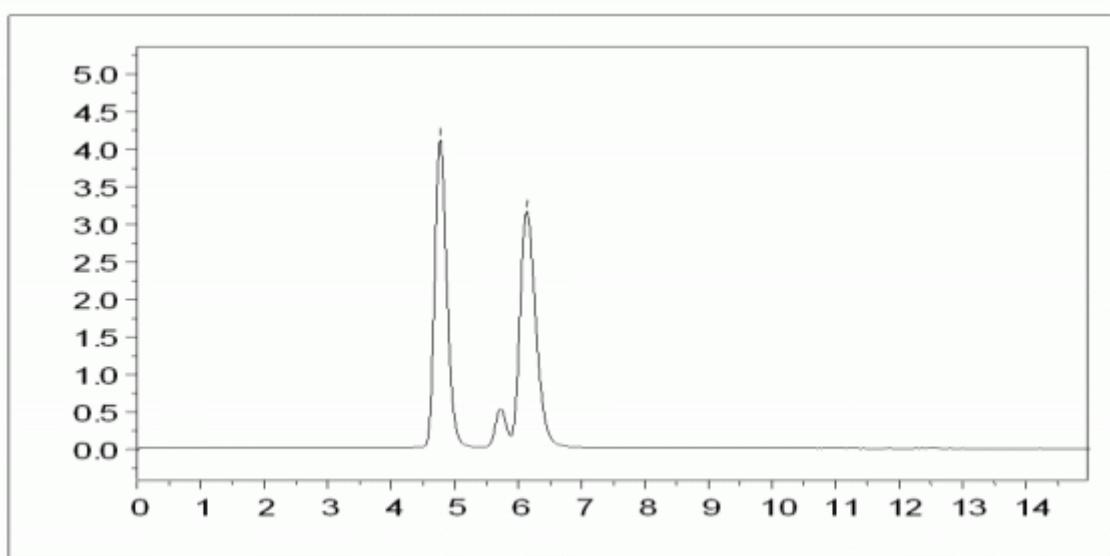

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® AD-RH		
Column Size (Particle Size)	4.6 x 150 mm (5 µm)		
Mobile Phase Composition	water / acetonitrile = 40 / 60		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	0.5	Temperature (°C)	25
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	4.8	0.33	1 - 2	3.40
2	6.1	0.71		