


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

alpha,beta-Isopropylidene glycol-gamma-phenyl ether

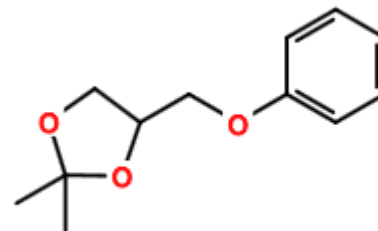

## CHROMATOGRAM

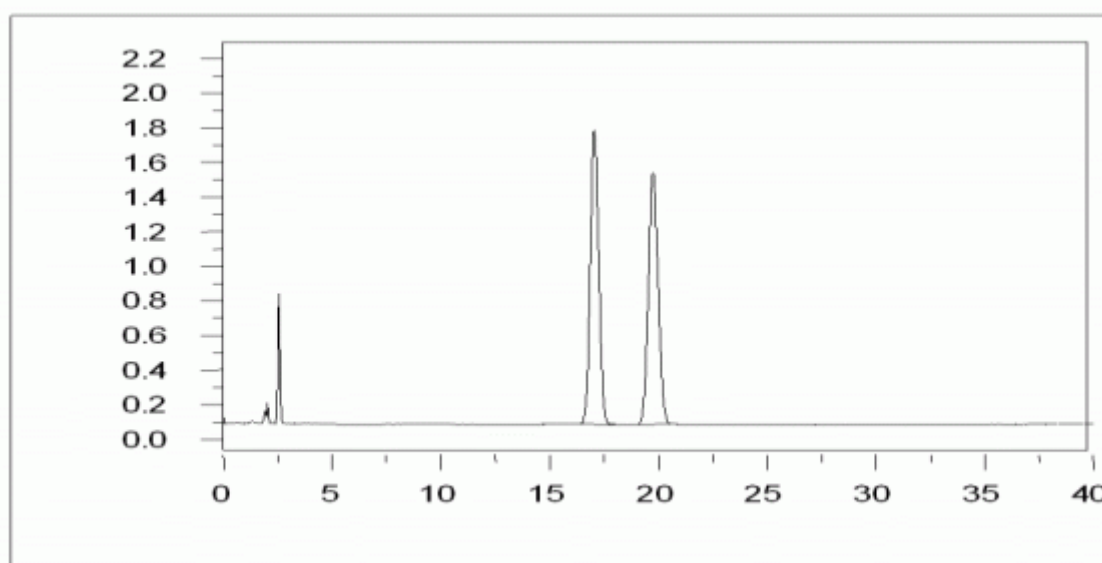

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® AS-RH		
Column Size (Particle Size)	4.6 x 150 mm (5 µm)		
Mobile Phase Composition	water / acetonitrile = 70 / 30		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	17.1	8.74	1 - 2	1.18
2	19.7	10.28		