


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Naproxen

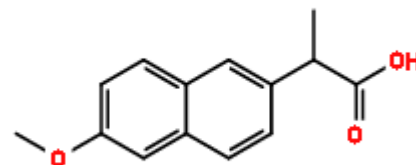

## CHROMATOGRAM

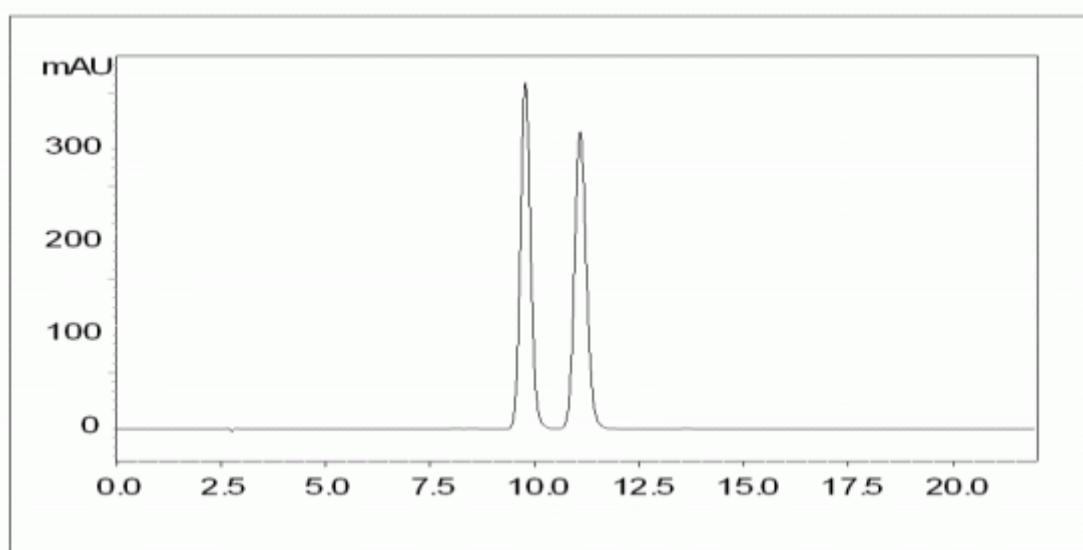

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® AS-RH		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	H3PO4 aq. pH2.0 / acetonitrile = 60 / 40		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	0.7	Temperature (°C)	25
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	9.8	2.81	1 - 2	2.51
2	11.1	3.33		