


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Perisoxal

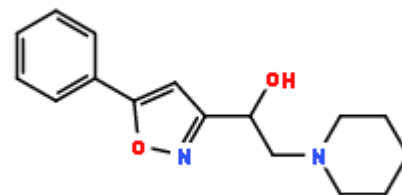

## CHROMATOGRAM

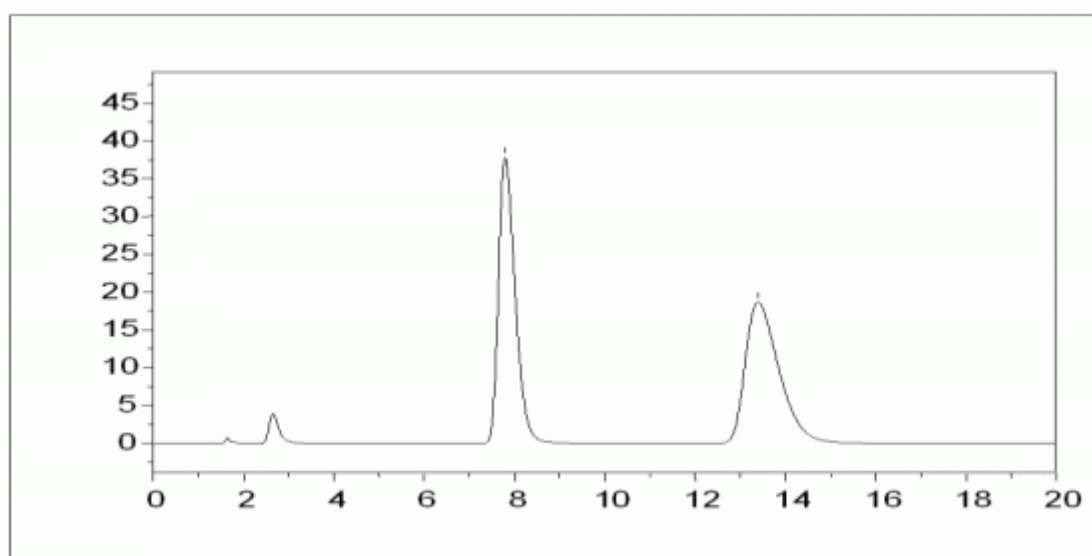

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® AD-RH		
Column Size (Particle Size)	4.6 x 150 mm (5 µm)		
Mobile Phase Composition	20mM borate buffer pH9.0 / acetonitrile = 30 / 70		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	7.8	5.19	1 - 2	5.51
2	13.4	9.63		