


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Pindolol

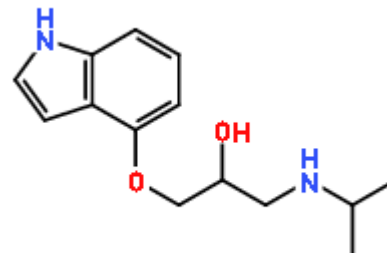

CHROMATOGRAM

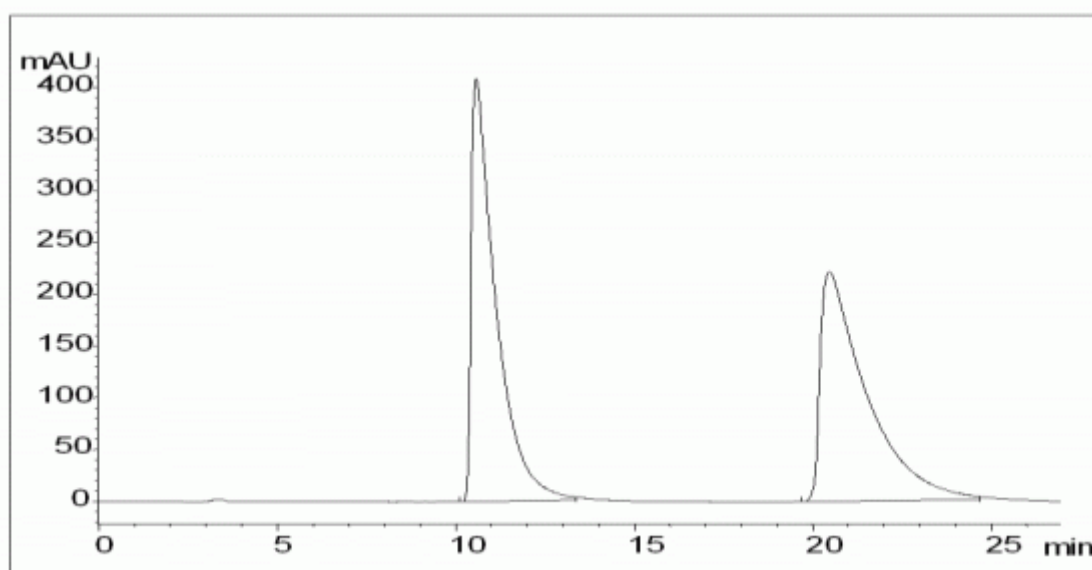

ANALYTICAL REPORT

Column Trade Name	CHIRALCEL® OD-RH		
Column Size (Particle Size)	4.6 x 150 mm (5 µm)		
Mobile Phase Composition	20mM borate buffer pH9.0 / acetonitrile = 60 / 40		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	0.5	Temperature (°C)	25
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 230 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	10.6	1.94	1 - 2	6.40
2	20.5	4.69		