


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Temazepam

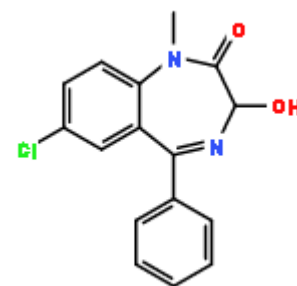

CHROMATOGRAM

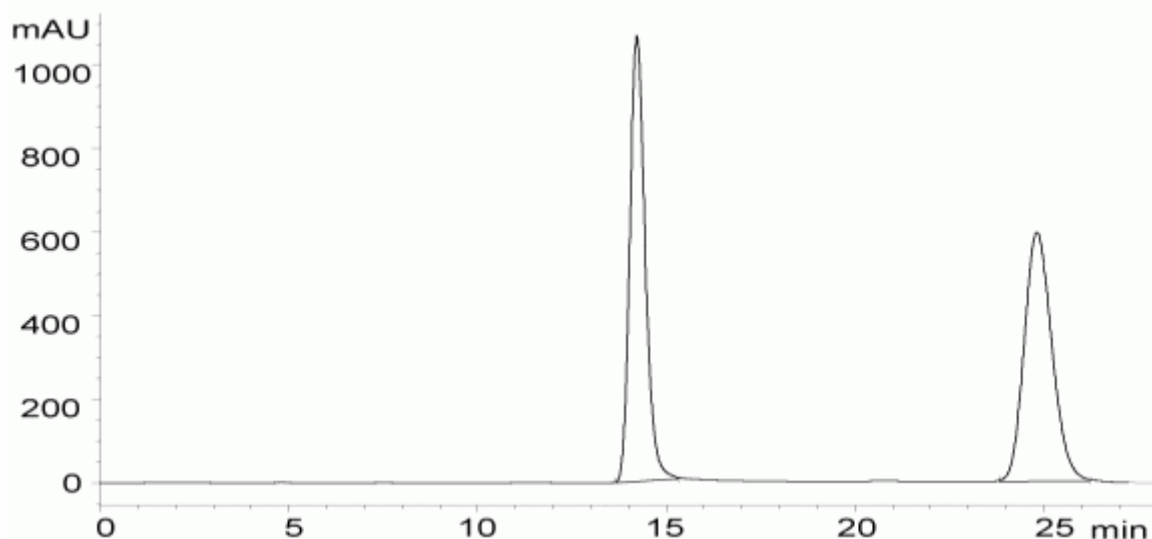

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® AS-RH		
Column Size (Particle Size)	4.6 x 150 mm (5 µm)		
Mobile Phase Composition	water / acetonitrile = 60 / 40		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	0.5	Temperature (°C)	25
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 250 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	14.2	2.55	1 - 2	10.00
2	24.8	5.20		