


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

1,1'-binaphthyl-2,2'-diyl hydrogenphosphate

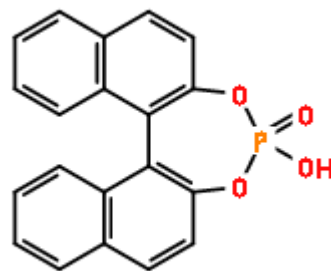

CHROMATOGRAM

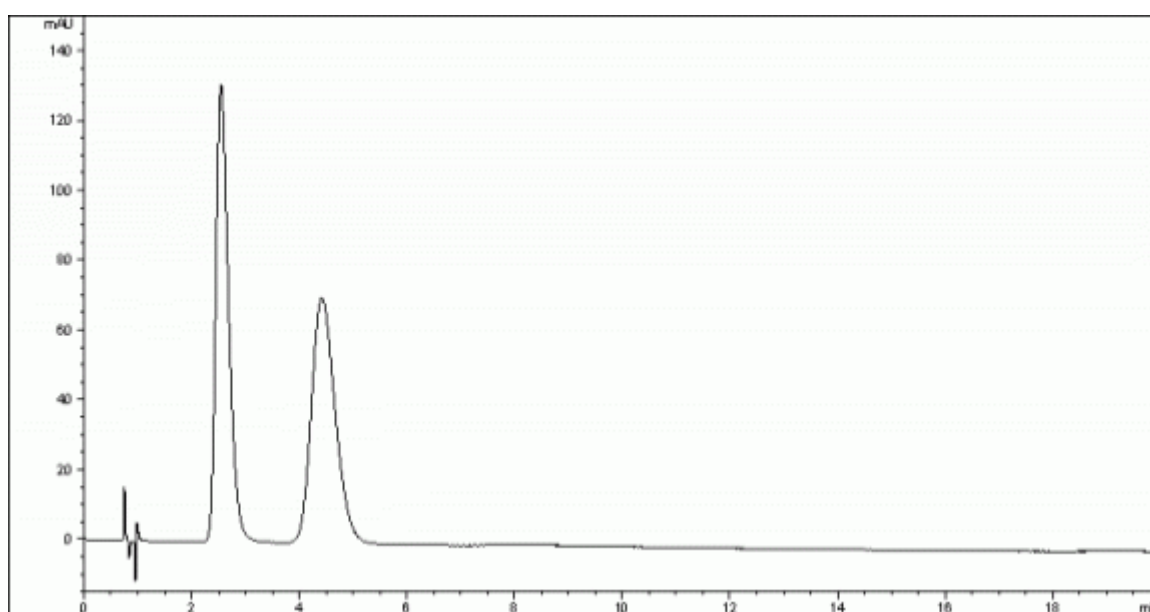

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® AGP		
Column Size (Particle Size)	4 x 100 mm (5 µm)		
Mobile Phase Composition	10mM NH4OAc aq. pH7.0 / acetonitrile = 85 / 15		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	0.9	Temperature (°C)	20
Injection Volume (µL)	0.0005 (mg)	Sample Concentration	
Detection	UV-VIS 225 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	2.6	1.51	1 - 2	2.22
2	4.4	3.36		