

APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

2-(2,4-Dichlorophenoxy)Propionic acid



CHROMATOGRAM

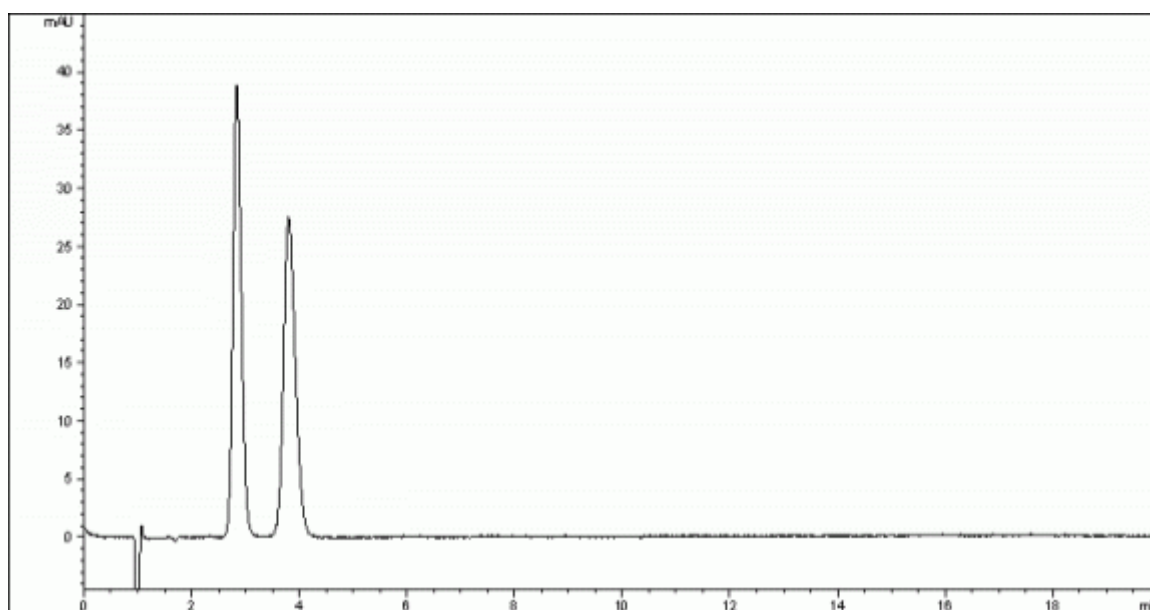

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® AGP		
Column Size (Particle Size)	4 x 100 mm (5 µm)		
Mobile Phase Composition	75mM NH4OAc aq. pH5.8 / 2-propanol = 95 / 5		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	0.9	Temperature (°C)	20
Injection Volume (µL)	0.0005 (mg)	Sample Concentration	
Detection	UV-VIS 225 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	2.9	1.80	1 - 2	2.82
2	3.8	2.76		