


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

alpha-Methoxyphenylacetic acid

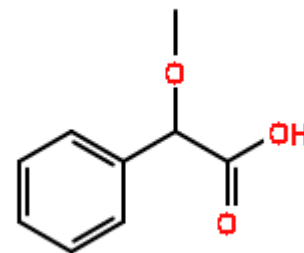

CHROMATOGRAM

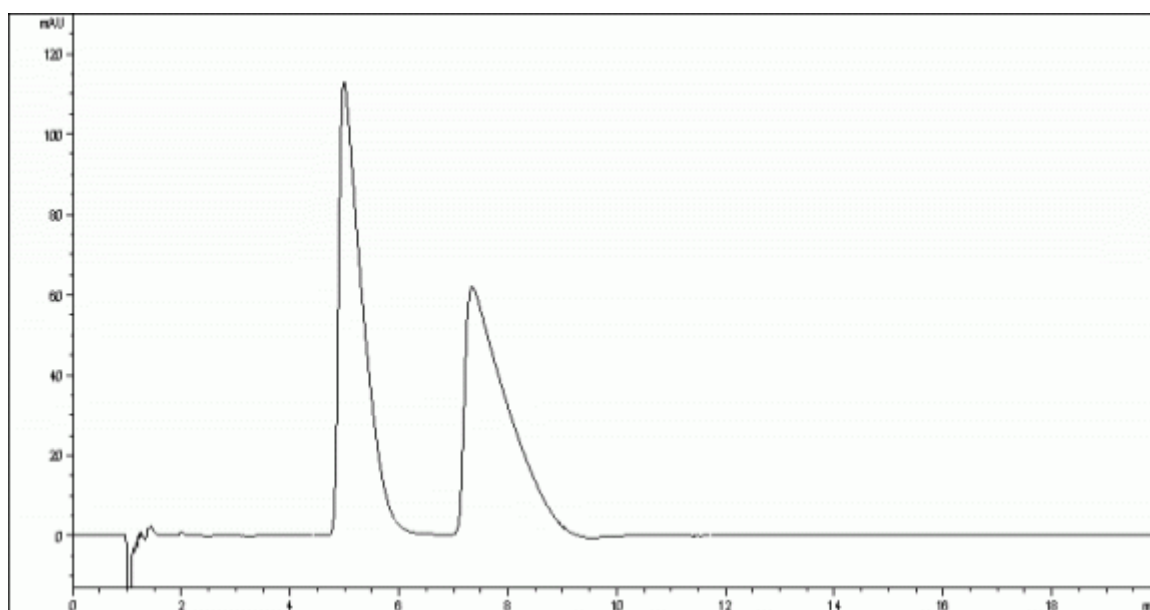

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® HSA		
Column Size (Particle Size)	4 x 100 mm (5 µm)		
Mobile Phase Composition	10mM NH ₄ OAc aq. pH7.0 / 2-propanol = 95 / 5		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	0.9	Temperature (°C)	20
Injection Volume (µL)	0.0025 (mg)	Sample Concentration	
Detection	UV-VIS 210 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	5.0	3.63	1 - 2	2.14
2	7.4	5.81		