


APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name

Benzoin

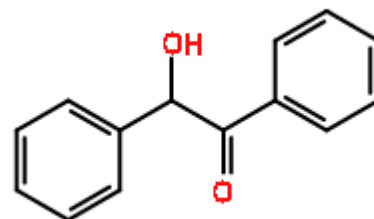

CHROMATOGRAM

NO CHROMATOGRAM

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® AGP		
Column Size (Particle Size)	4 x 100 mm (5 µm)		
Mobile Phase Composition	10mM NH ₄ OAc aq. pH7.0 / 2-propanol = 95 / 5		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	0.9	Temperature (°C)	20
Injection Volume (µL)	0.0005 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	5.1	3.79	1 - 2	1.26
2	6.1	4.77		1.96