


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Etodolac

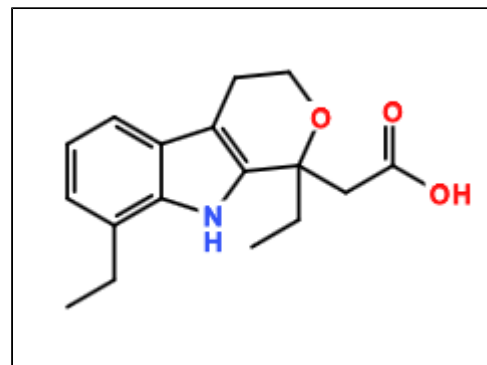

CHROMATOGRAM

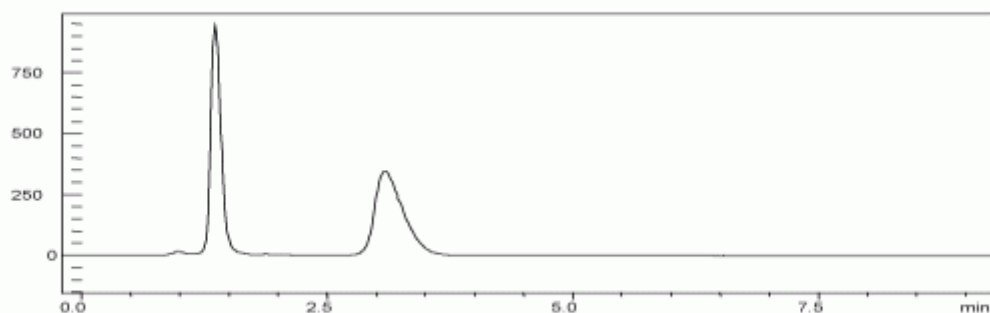

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® AGP		
Column Size (Particle Size)	4 x 100 mm (5 µm)		
Mobile Phase Composition	5mM NH ₄ OAc aq. pH7.0 / acetonitrile = 85 / 15		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	0.8	Temperature (°C)	25
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 220 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	1.3	0.40	1 - 2	3.52
2	3.2	2.57		