


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

O-Methoxymandelic acid

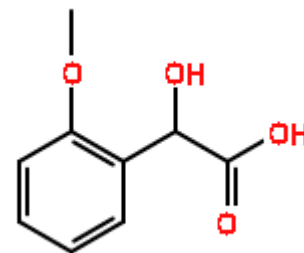

## CHROMATOGRAM

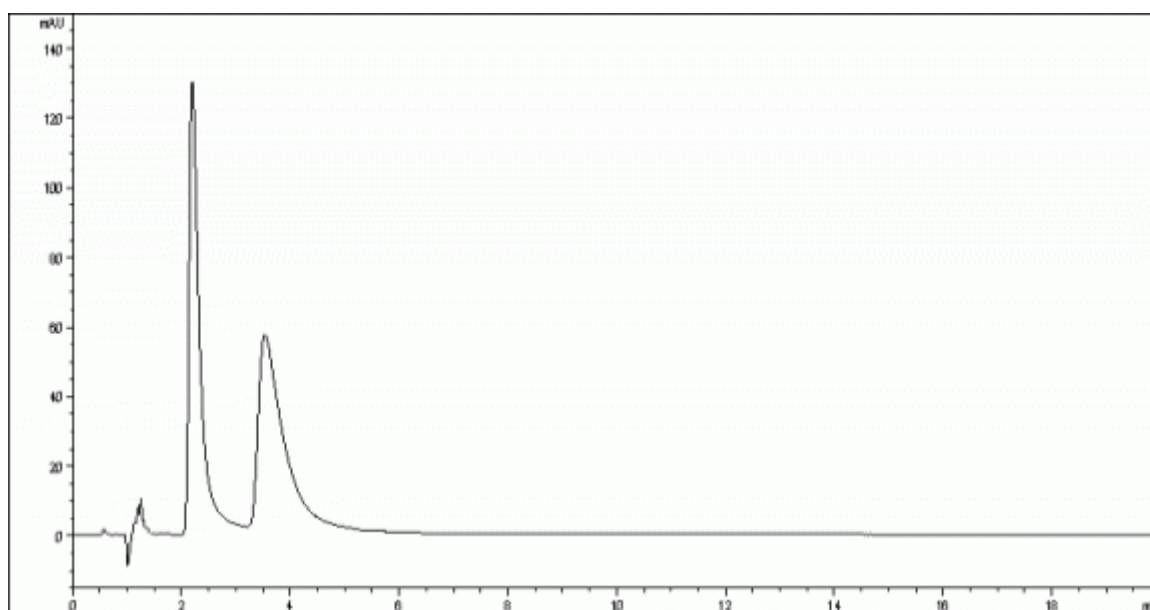

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK <sup>®</sup> HSA		
Column Size (Particle Size)	4 x 100 mm (5 μm)		
Mobile Phase Composition	10mM NH <sub>4</sub> OAc aq. pH7.0 / 2-propanol = 85 / 15		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	0.9	Temperature (°C)	20
Injection Volume (μL)	0.0025 (mg)	Sample Concentration	
Detection	UV-VIS 225 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	2.2	1.04	1 - 2	2.45
2	3.5	2.28		