


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

3-oxo-1-Indancarboxylic acid

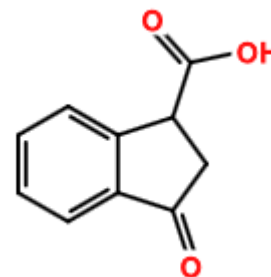

## CHROMATOGRAM

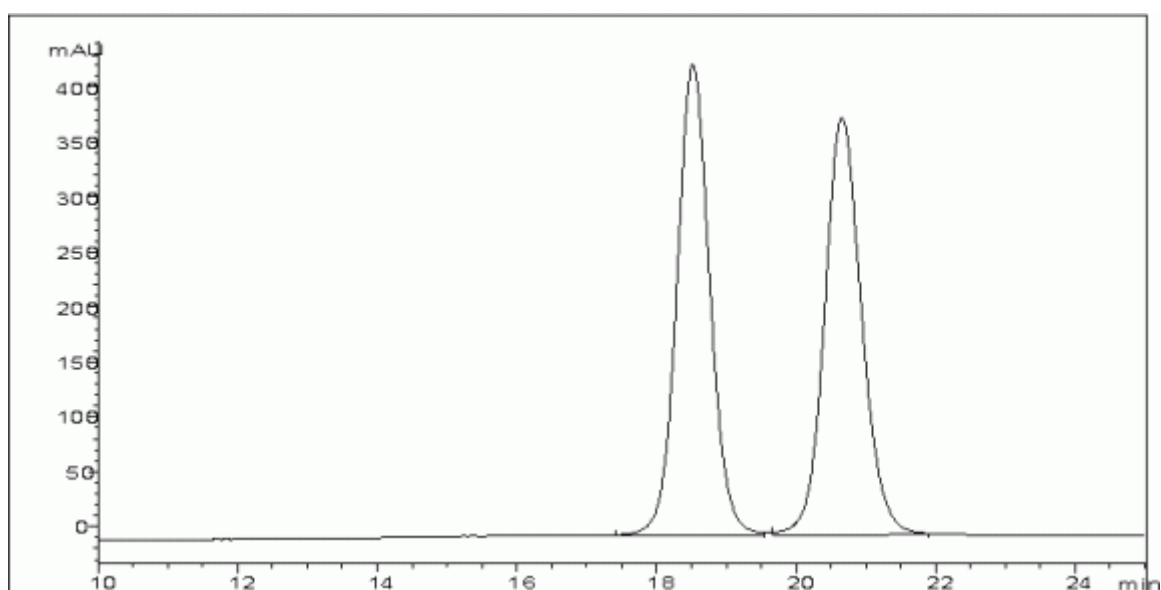

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® QN-AX		
Column Size (Particle Size)	4.6 x 150 mm (5 µm)		
Mobile Phase Composition	methanol / 0.2mM acetic acid aq. pH6.0 = 80 / 20		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	18.5	8.96	1 - 2	2.41
2	20.7	10.11		