



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	3-(1-Naphthyl)-DL-alanine
---------------	---------------------------

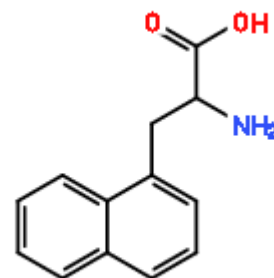

CHROMATOGRAM

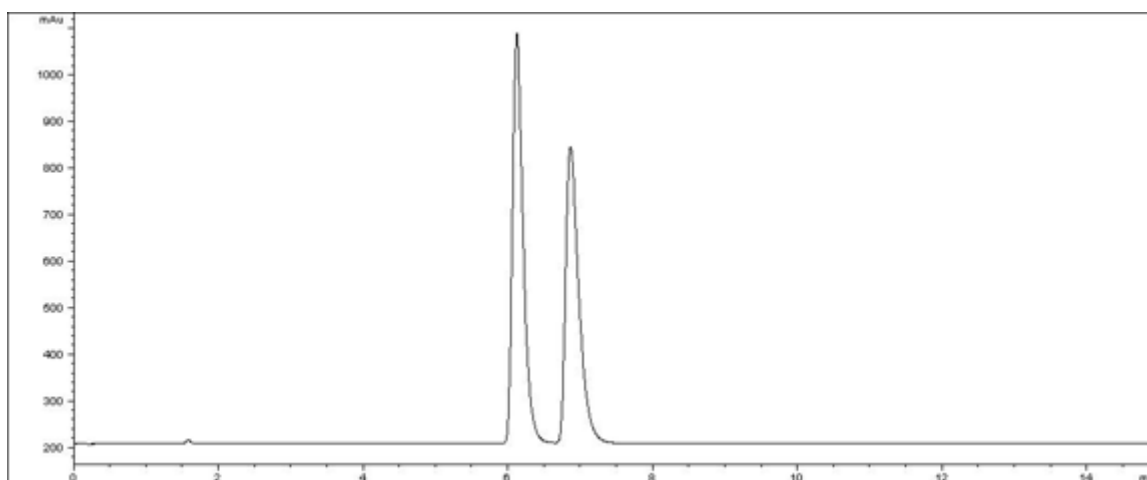

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® ZWIX(+)
Column Size / Particle Size	150 x 3mm / 3µm
Mobile Phase Composition	50mM formic acid + 25mM diethylamine in methanol / acetonitrile / water = 49 / 49 / 2
Chromatographic Mode	Polar Phase

Flow rate (ml/min)	0.5	Injection amount (mg)	0.005
Temperature (°C)	10	Pressure (MPa)	
Detection	ELSD	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	6.12	2.36	1 - 2	1.17
2	6.86	2.77		2.56