



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	3-(2-Thienyl)-DL-alanine
---------------	--------------------------

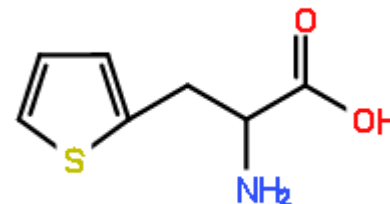

CHROMATOGRAM

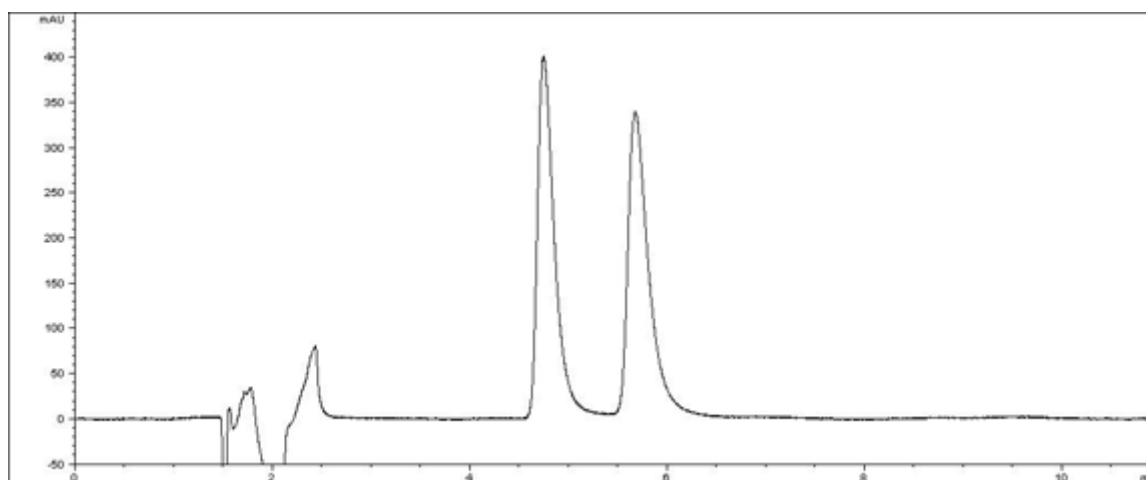

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® ZWIX(-) 150 x 3mm / 3µm
Mobile Phase Composition	50mM formic acid + 25mM diethylamine in methanol / tetrahydrofuran / water = 49 / 49 / 2
Chromatographic Mode	Polar Phase

Flow rate (ml/min)	0.5	Injection amount (mg)	0.0025
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 230 nm	Elution order	D/L

No.	Rt	K'	Alpha Value	Resolution
1	4.75	1.93	1 - 2	1.29
2	5.67	2.5		2.8