


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Metalaxyl

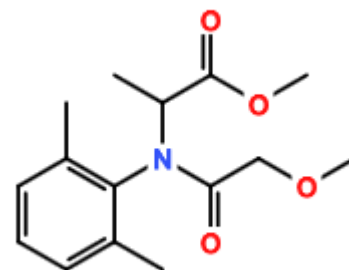

CHROMATOGRAM

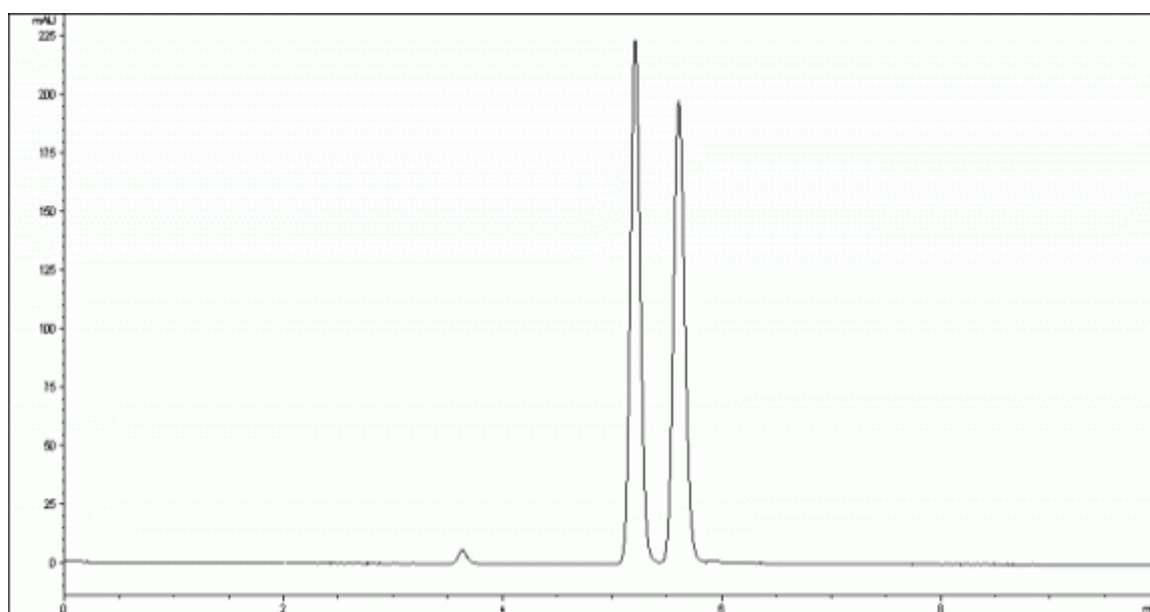

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALCEL® OZ-3 150 x 4.6mm / 3µm
Mobile Phase Composition	ethanol = 100
Chromatographic Mode	Polar Phase

Flow rate (ml/min)	0.5	Injection amount (mg)	0.005
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 230 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	5.21	2.08	1 - 2 1.11	2.26
2	5.61	2.32		