



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	Sertraline tetralone
---------------	----------------------

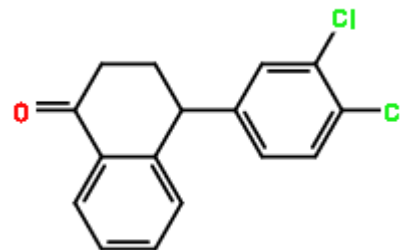

CHROMATOGRAM

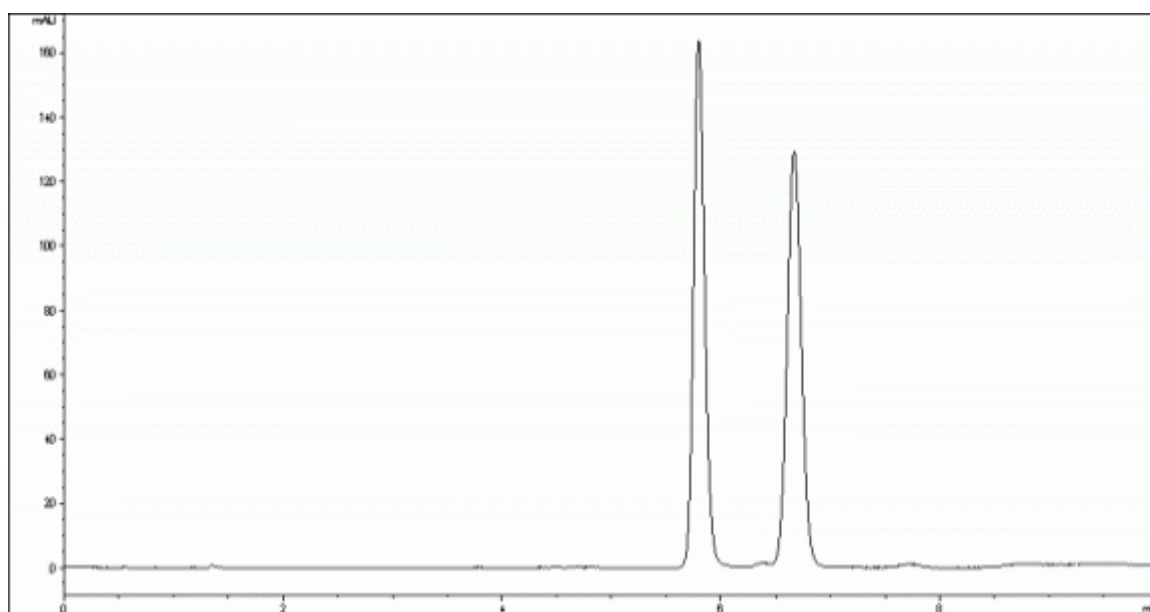

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® AY-3 150 x 4.6mm / 3µm
Mobile Phase Composition	ethanol = 100
Chromatographic Mode	Polar Phase

Flow rate (ml/min)	0.5	Injection amount (mg)	0.005
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 270 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	5.8	2.28	1 - 2	4.15
2	6.67	2.77		