


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Thalidomide

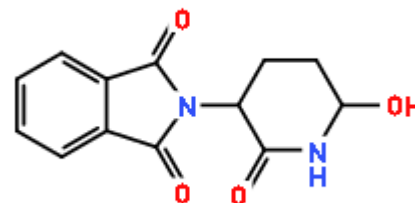

## CHROMATOGRAM

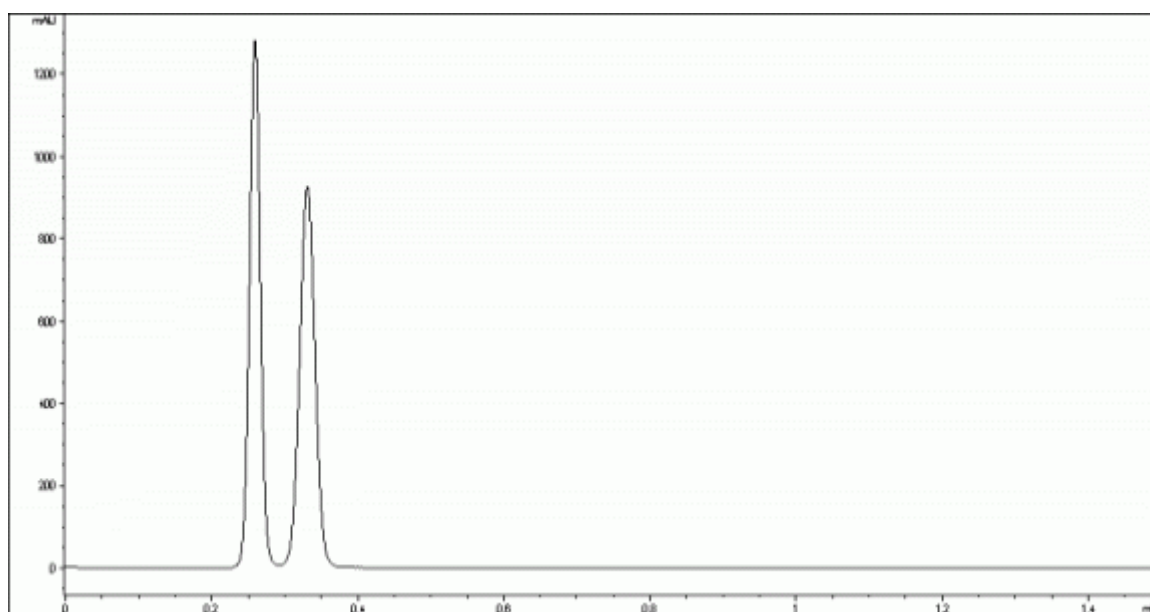

## ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALCEL® OZ-3 50 x 4.6mm / 3µm
Mobile Phase Composition	acetonitrile = 100
Chromatographic Mode	Polar Phase

Flow rate (ml/min)	4	Injection amount (mg)	0.005
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 230 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	0.26	0.69	1 - 2	2.19
2	0.33	1.16		