



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	Wieland Miescher Ketone
---------------	-------------------------

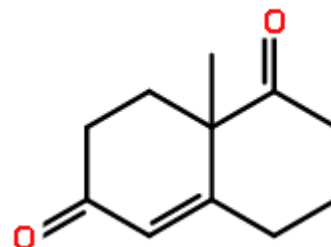

CHROMATOGRAM

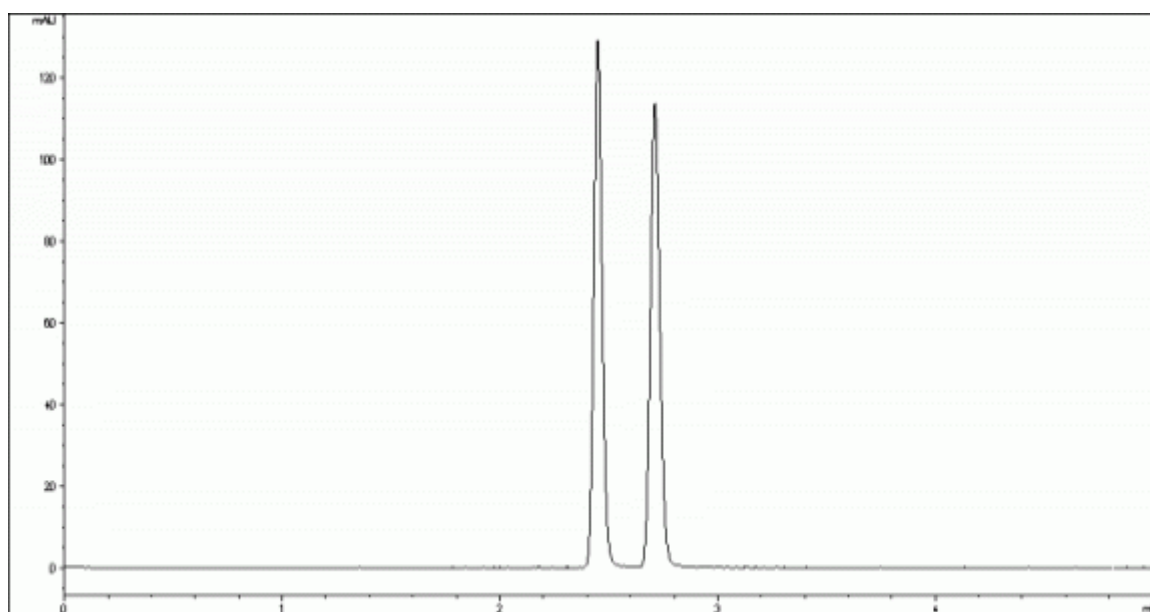

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALCEL® OZ-3 150 x 4.6mm / 3µm
Mobile Phase Composition	acetonitrile = 100
Chromatographic Mode	Polar Phase

Flow rate (ml/min)	1	Injection amount (mg)	0.005
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 270 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	2.45	0.45	1 - 2 1.34	3.54
2	2.71	0.6		