



## APPLICATION RESULTS

## COMPOUND IDENTIFICATION

Compound Name	5-(p-Methylphenyl)-5-phenylhydantoin
---------------	--------------------------------------

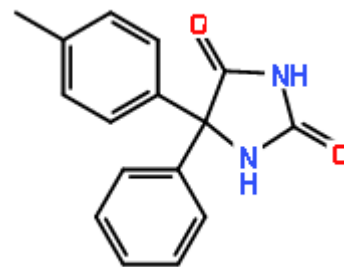

## CHROMATOGRAM

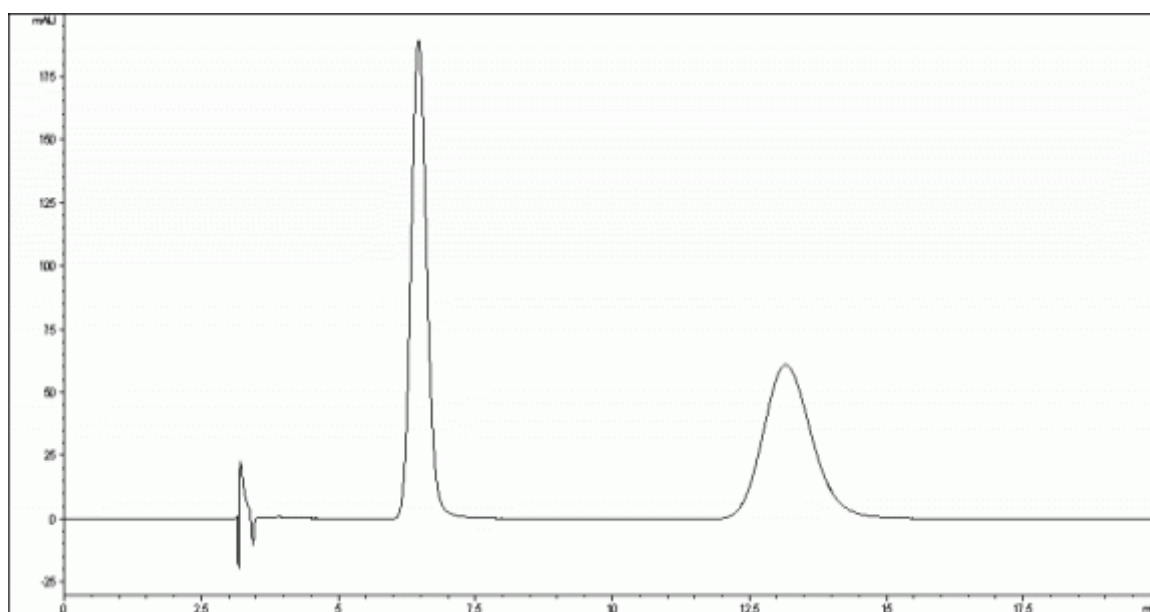

## ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® AZ-H 250 x 4.6mm / 5µm
Mobile Phase Composition	acetonitrile / 2-propanol = 95 / 5
Chromatographic Mode	Polar Phase

Flow rate (ml/min)	1	Injection amount (mg)	0.005
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 230 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	6.47	1.18	1 - 2 2.91	6.28
2	13.16	3.43		