



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	Nomifensine
---------------	-------------

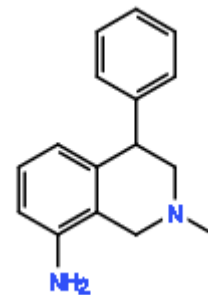

CHROMATOGRAM

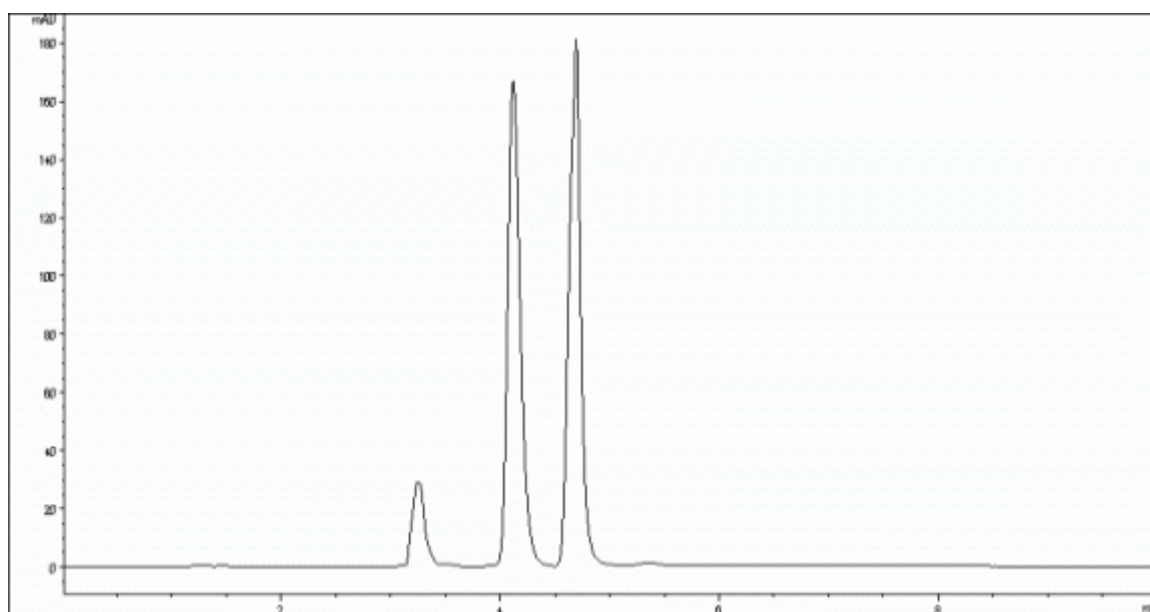

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALCEL® OZ-H 250 x 4.6mm / 5µm
Mobile Phase Composition	acetonitrile / diethylamine = 100 / 0.1
Chromatographic Mode	Polar Phase

Flow rate (ml/min)	1	Injection amount (mg)	0.01
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 270 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	4.12	0.37	1 - 2	1.52
2	4.7	0.56		2.9