


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Ornidazole

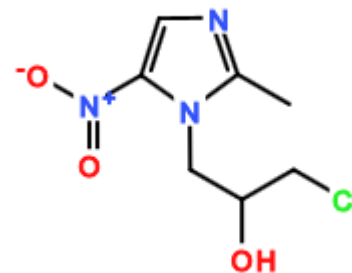

## CHROMATOGRAM

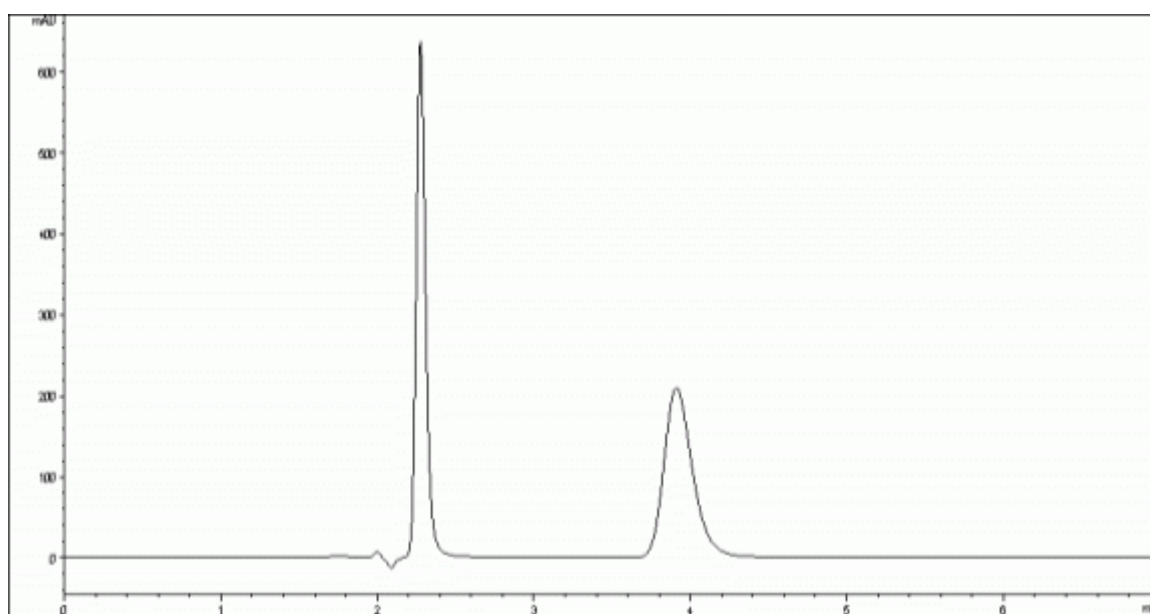

## ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® AY-H 150 x 4.6mm / 5µm
Mobile Phase Composition	ethanol / diethylamine = 100 / 0.1
Chromatographic Mode	Polar Phase

Flow rate (ml/min)	1	Injection amount (mg)	0.005
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 230 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	2.28	0.22	1 - 2	4.94
2	3.91	1.1		7.66