



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	Sertraline tetralone
---------------	----------------------

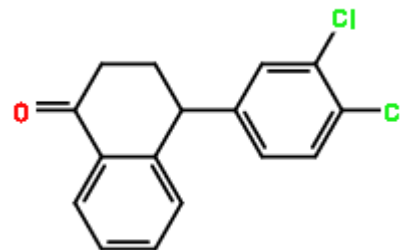

CHROMATOGRAM

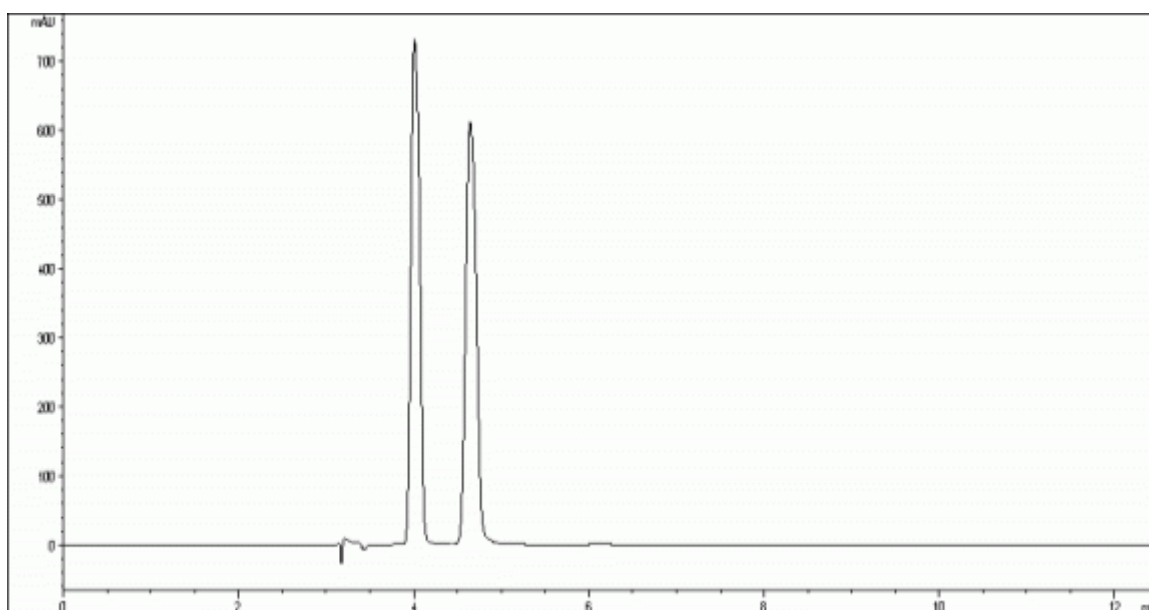

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® AZ-H 250 x 4.6mm / 5µm
Mobile Phase Composition	acetonitrile / 2-propanol = 95 / 5
Chromatographic Mode	Polar Phase

Flow rate (ml/min)	1	Injection amount (mg)	0.005
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 254 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	4.01	0.35	1 - 2	3.02
2	4.65	0.57		