


APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	Aminogluthethimide
---------------	--------------------



CHROMATOGRAM

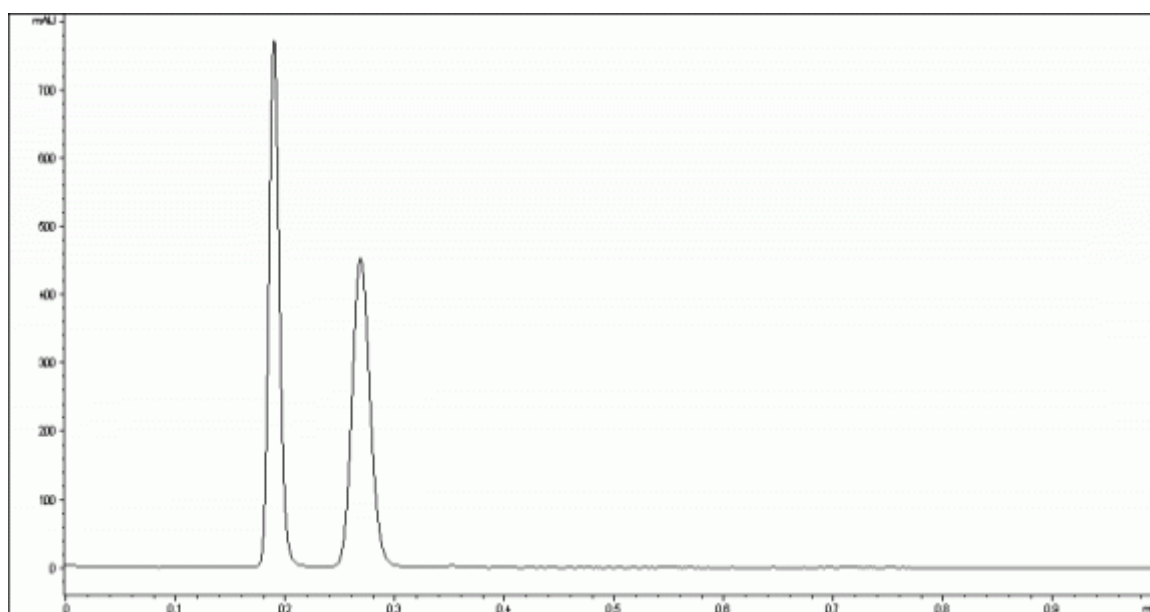

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALCEL® OJ-3 50 x 4.6mm / 3µm
Mobile Phase Composition	methanol = 100
Chromatographic Mode	Polar Phase

Flow rate (ml/min)	5	Injection amount (mg)	0.004
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 254 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	0.19	0.59	1 - 2 2.12	3.24
2	0.27	1.25		