



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	N-alpha-Benzoyl-DL-arginine-beta-naphtylamide
---------------	---

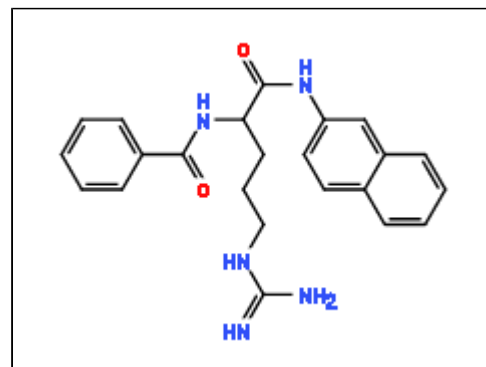

CHROMATOGRAM

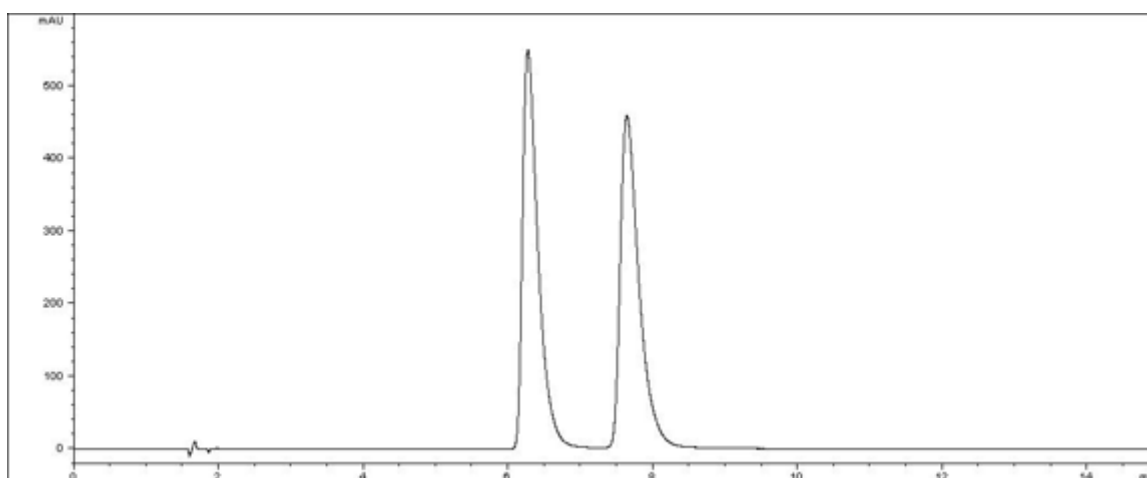

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® ZWIX(+) 150 x 3mm / 3µm
Mobile Phase Composition	50mM formic acid + 25mM diethylamine in methanol / water = 90 / 10
Chromatographic Mode	Polar Phase

Flow rate (ml/min)	0.5	Injection amount (mg)	0.0025
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 254 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	6.27	2.79	1 - 2	3.31
2	7.64	3.62		