



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	Wieland Miescher Ketone
---------------	-------------------------

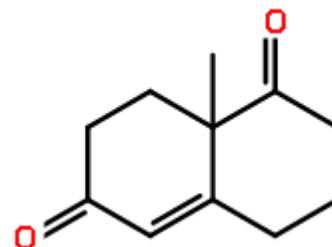

CHROMATOGRAM

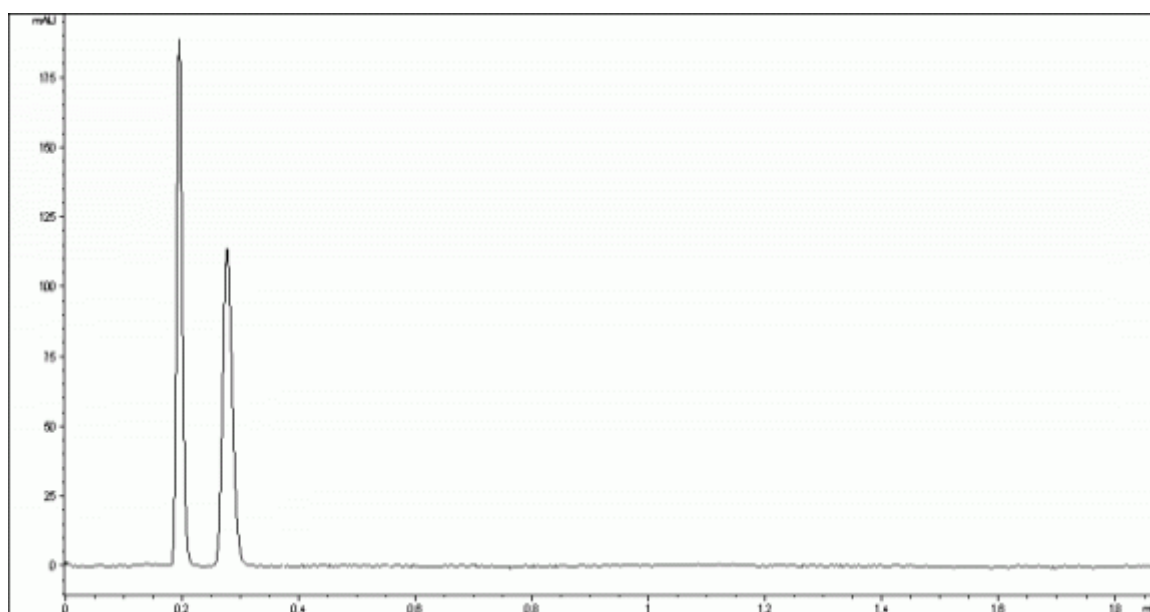

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® AS-3 50 x 4.6mm / 3µm
Mobile Phase Composition	acetonitrile = 100
Chromatographic Mode	Polar Phase

Flow rate (ml/min)	5	Injection amount (mg)	0.01
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 270 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	0.2	0.75	1 - 2	3.35
2	0.28	1.48		