



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	(±)-2,3-O-Benzylidene threitol
---------------	--------------------------------

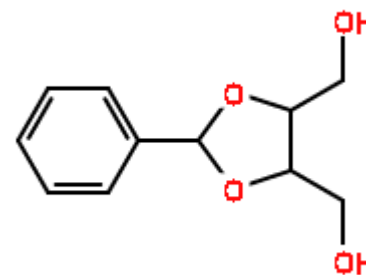

CHROMATOGRAM

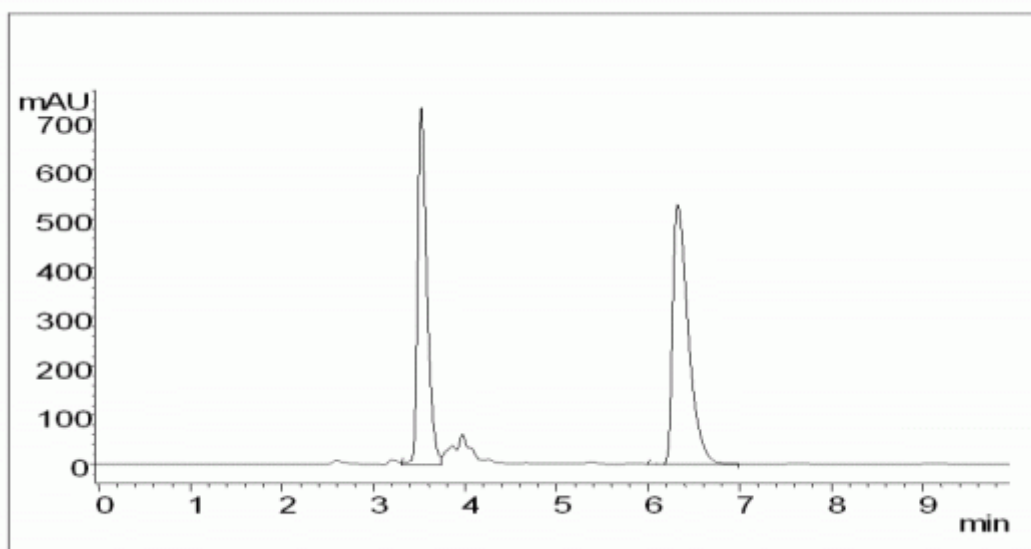

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® AD-H 250 x 4.6mm / 5µm
Mobile Phase Composition	methanol = 100
Chromatographic Mode	Polar Phase

Flow rate (ml/min)	1	Injection amount (mg)	0.005
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 220 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	3.52	0.17	1 - 2	11.83
2	6.32	1.11		