



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	N-benzoyl-DL-phenylalanine
---------------	----------------------------

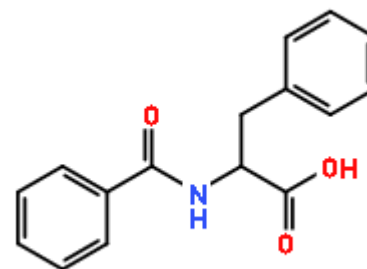

CHROMATOGRAM

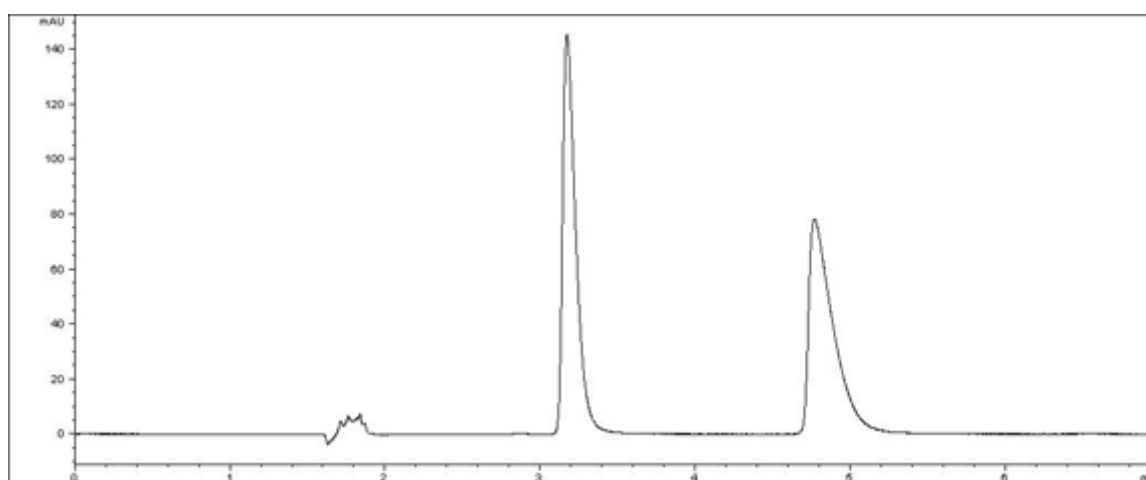

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® ZWIX(+)
Column Size / Particle Size	150 x 3mm / 3µm
Mobile Phase Composition	25mM acetic acid in methanol / water = 98 / 2
Chromatographic Mode	Polar Phase

Flow rate (ml/min)	0.5	Injection amount (mg)	0.0025
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 254 nm	Elution order	D/L

No.	Rt	K'	Alpha Value	Resolution
1	3.17	0.93	1 - 2	7.54
2	4.77	1.9		