



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	N-CBZ-DL-asparagine
---------------	---------------------

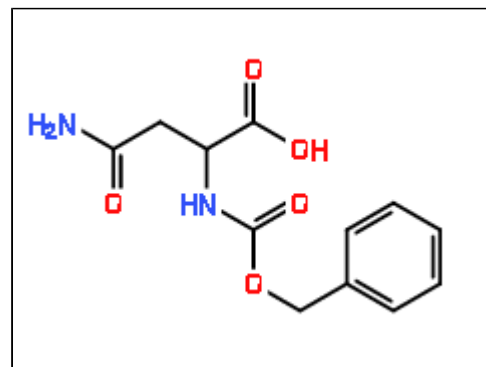

CHROMATOGRAM

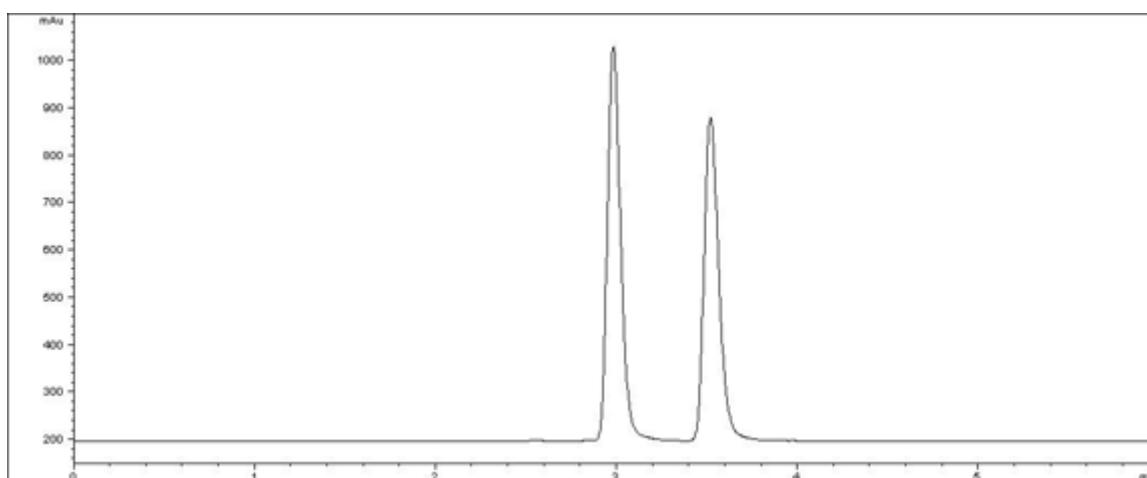

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® ZWIX(+) 150 x 3mm / 3µm
Mobile Phase Composition	50mM formic acid + 25mM diethylamine in methanol / acetonitrile / water = 49 / 49 / 2
Chromatographic Mode	Polar Phase

Flow rate (ml/min)	0.5	Injection amount (mg)	0.0025
Temperature (°C)	25	Pressure (MPa)	
Detection	ELSD	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	2.98	0.7	1 - 2	3.88
2	3.52	1.01		