



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	Chlormezanone
---------------	---------------

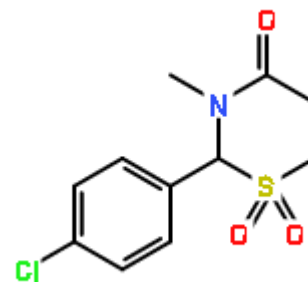

CHROMATOGRAM

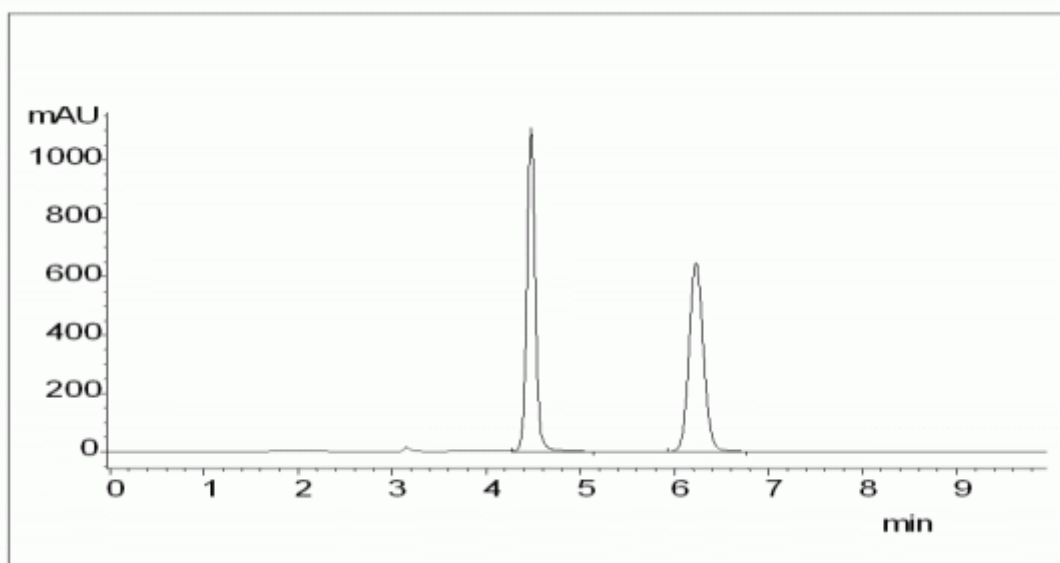

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® AS-H 250 x 4.6mm / 5µm
Mobile Phase Composition	methanol = 100
Chromatographic Mode	Polar Phase

Flow rate (ml/min)	1	Injection amount (mg)	0.005
Temperature (°C)	30	Pressure (MPa)	
Detection	UV 210 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	4.48	0.49	1 - 2 2.19	8.2
2	6.23	1.08		