



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	N-CBZ-DL-glutamine
---------------	--------------------

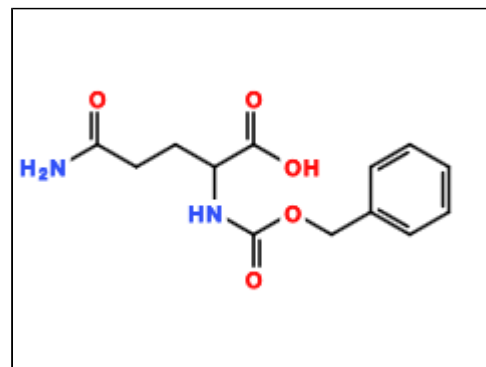

CHROMATOGRAM

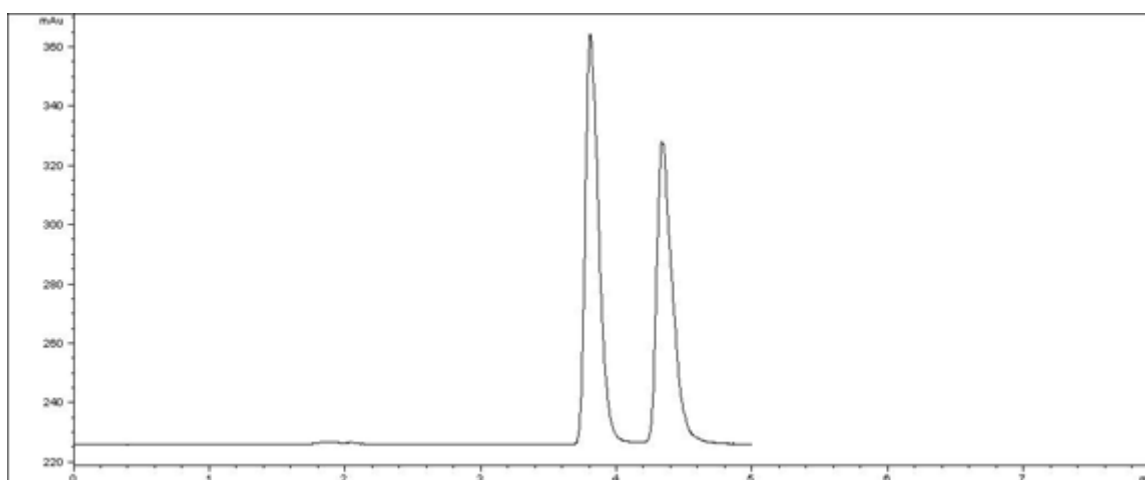

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® ZWIX(+) 150 x 3mm / 3µm
Mobile Phase Composition	25mM acetic acid in methanol / water = 98 / 2
Chromatographic Mode	Polar Phase

Flow rate (ml/min)	0.5	Injection amount (mg)	0.0025
Temperature (°C)	10	Pressure (MPa)	
Detection	ELSD	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	3.81	1.08	1 - 2	2.71
2	4.33	1.37		