



## APPLICATION RESULTS

## COMPOUND IDENTIFICATION

Compound Name	Mephobarbital
---------------	---------------

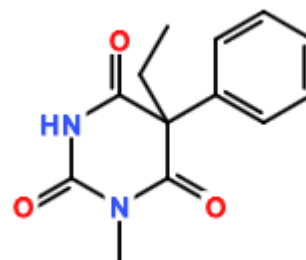

## CHROMATOGRAM

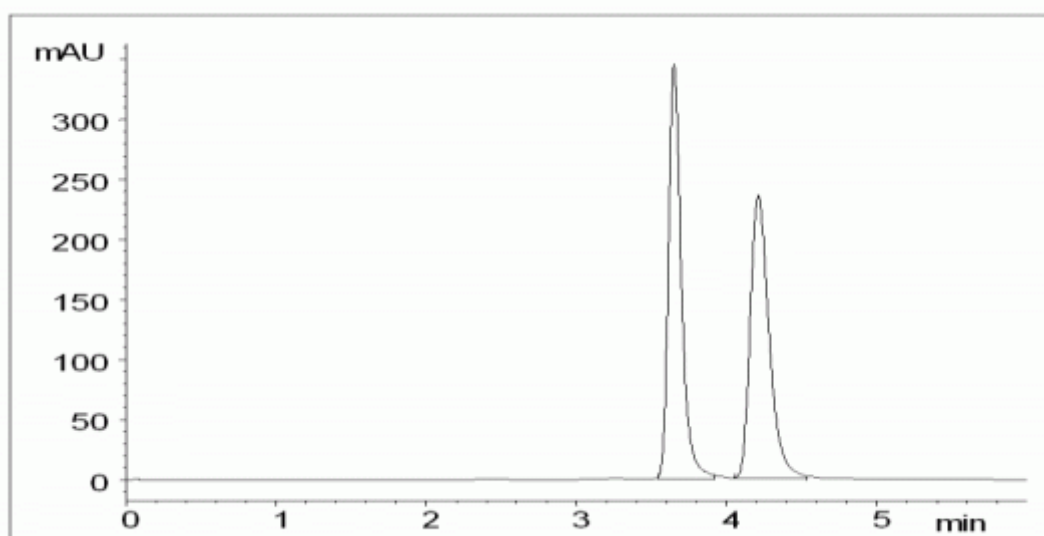

## ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® AS-H 250 x 4.6mm / 5µm
Mobile Phase Composition	methanol = 100
Chromatographic Mode	Polar Phase

Flow rate (ml/min)	1	Injection amount (mg)	0.01
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 210 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	3.65	0.22	1 - 2	2.9
2	4.21	0.4		