



## APPLICATION RESULTS

## COMPOUND IDENTIFICATION

Compound Name	N-CBZ-DL-Norvaline
---------------	--------------------

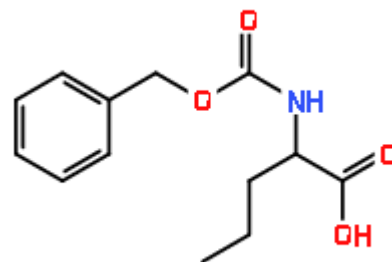

## CHROMATOGRAM

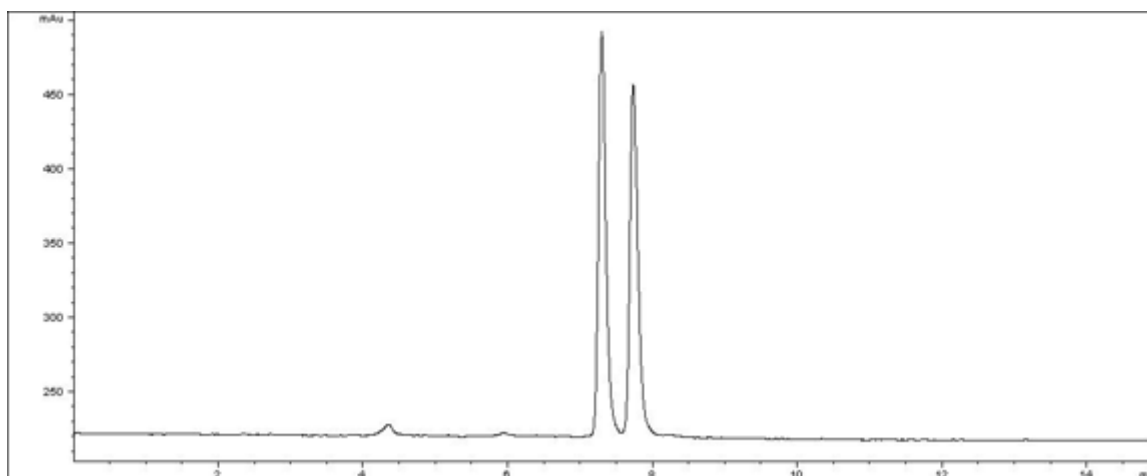

## ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® ZWIX(+) 250 x 4mm / 3µm
Mobile Phase Composition	25mM acetic acid in methanol / water = 98 / 2
Chromatographic Mode	Polar Phase

Flow rate (ml/min)	0.5	Injection amount (mg)	0.0015
Temperature (°C)	10	Pressure (MPa)	
Detection	ELSD	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	7.3	1.57	1 - 2	2.26
2	7.73	1.72		