



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	N-CBZ-DL-Valine
---------------	-----------------

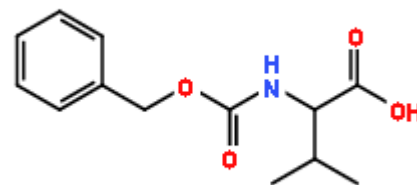

CHROMATOGRAM

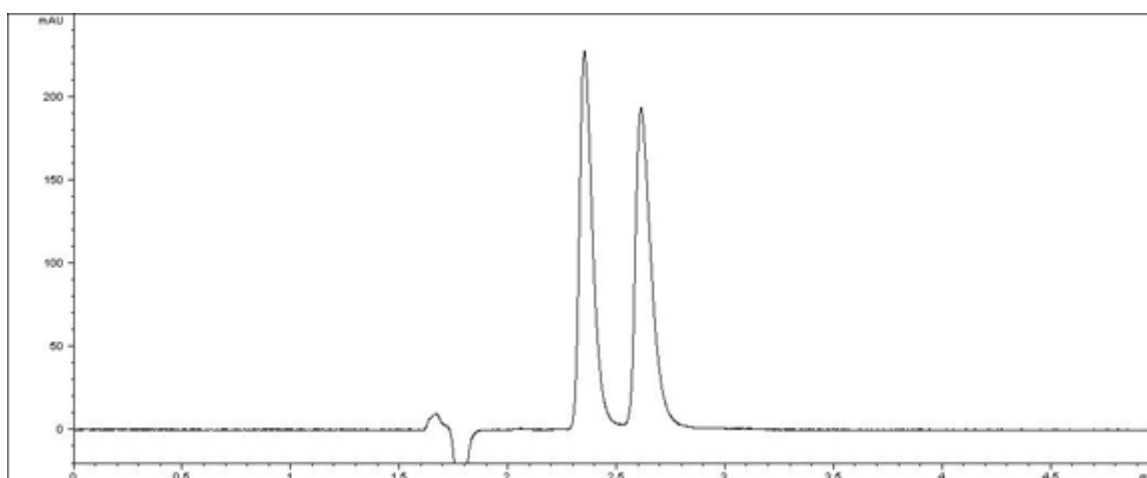

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® ZWIX(+) 150 x 3mm / 3µm
Mobile Phase Composition	25mM acetic acid in methanol / water = 98 / 2
Chromatographic Mode	Polar Phase

Flow rate (ml/min)	0.5	Injection amount (mg)	0.0025
Temperature (°C)	10	Pressure (MPa)	
Detection	UV 218 nm	Elution order	D/L

No.	Rt	K'	Alpha Value	Resolution
1	2.35	0.41	1 - 2	2.19
2	2.61	0.56		