



## APPLICATION RESULTS

## COMPOUND IDENTIFICATION

Compound Name	N-FMOC-DL-Phenylalanine
---------------	-------------------------

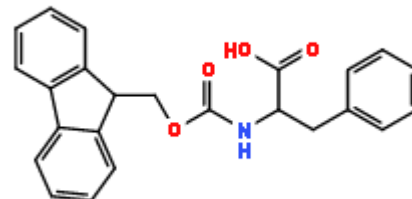

## CHROMATOGRAM

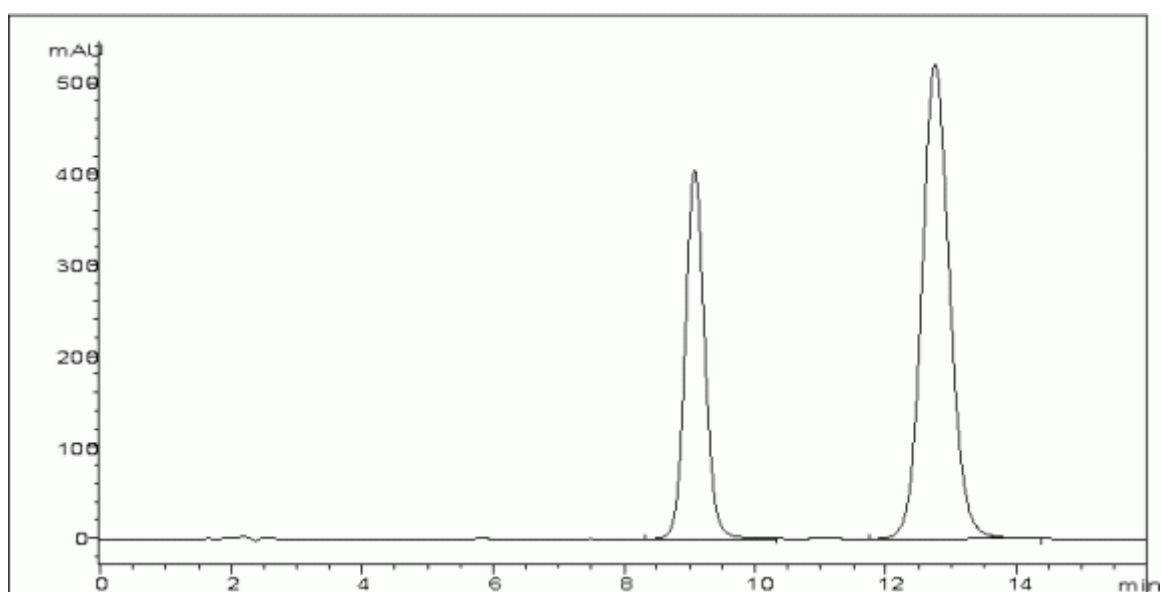

## ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® QN-AX 150 x 4.6mm / 5µm
Mobile Phase Composition	methanol / acetic acid / NH <sub>4</sub> OAc = 98 / 2 / 0.5
Chromatographic Mode	Polar Phase

Flow rate (ml/min)	1	Injection amount (mg)	0.01
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 254 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	9.08	3.88	1 - 2	5.71
2	12.75	5.86		