



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	(±)-2,3-O-Benzylidene threitol
---------------	--------------------------------

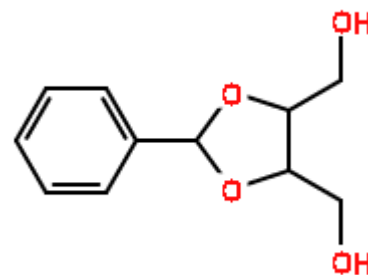

CHROMATOGRAM

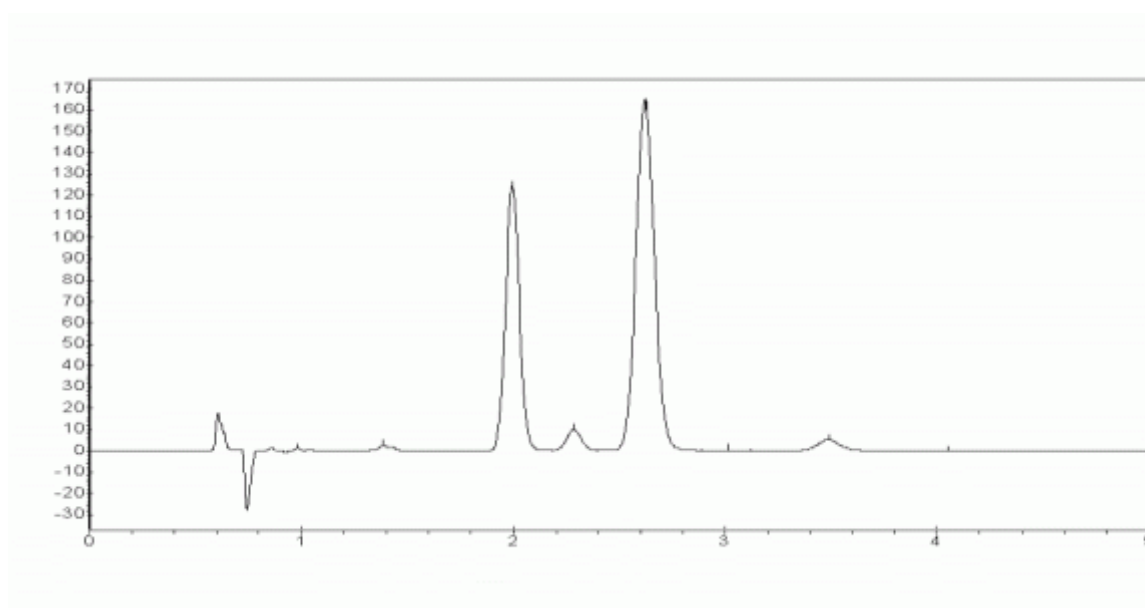

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® IC 250 x 4.6mm / 5µm
Mobile Phase Composition	CO ₂ / 2-propanol = 80 / 20
Chromatographic Mode	SFC

Flow rate (ml/min)	5	Injection amount (mg)	0.01
Temperature (°C)	35	Pressure (MPa)	10
Detection	UV 210 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	2	2.33	1 - 2	1.44
2	2.62	3.37		