



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	(4R,5S/4S,5R)-1,5-Dimethyl-4-phenyl-2-imidazolidinone
---------------	---

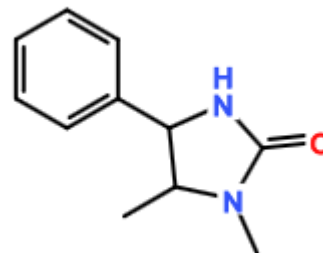

CHROMATOGRAM

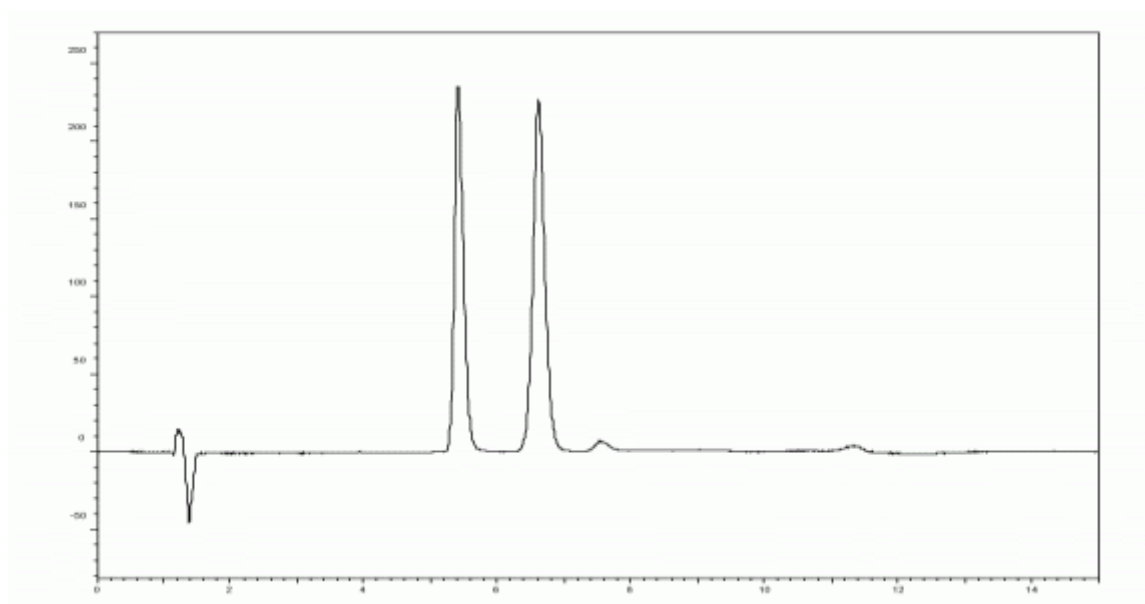

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® AY-H 250 x 4.6mm / 5µm
Mobile Phase Composition	CO ₂ / 2-propanol = 85 / 15
Chromatographic Mode	SFC

Flow rate (ml/min)	3	Injection amount (mg)	
Temperature (°C)	RT	Pressure (MPa)	12
Detection	UV 210 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	5.41	4.41	1 - 2	4.71
2	6.62	5.62		