



## APPLICATION RESULTS

## COMPOUND IDENTIFICATION

Compound Name	1-(Benzyloxycarbonyl)-2-tert-butyl-3-methyl-4-imidazolidinone
---------------	---

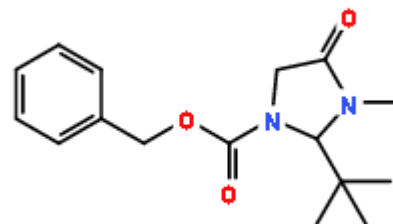

## CHROMATOGRAM

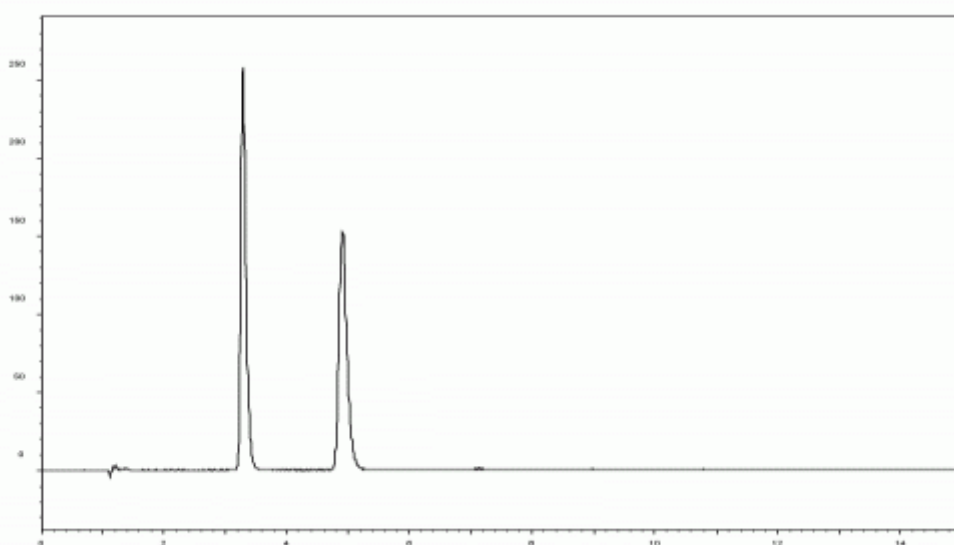

## ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALCEL® OZ-H 250 x 4.6mm / 5µm
Mobile Phase Composition	CO <sub>2</sub> / methanol = 85 / 15
Chromatographic Mode	SFC

Flow rate (ml/min)	3	Injection amount (mg)	
Temperature (°C)	RT	Pressure (MPa)	12
Detection	UV 210 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	3.29	2.29	1 - 2	8.37
2	4.91	3.91		