



## APPLICATION RESULTS

## COMPOUND IDENTIFICATION

Compound Name	Benzoin ethyl ether
---------------	---------------------

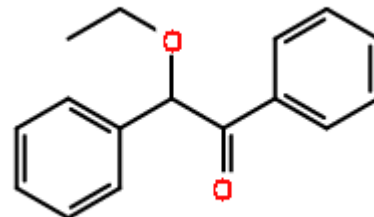

## CHROMATOGRAM

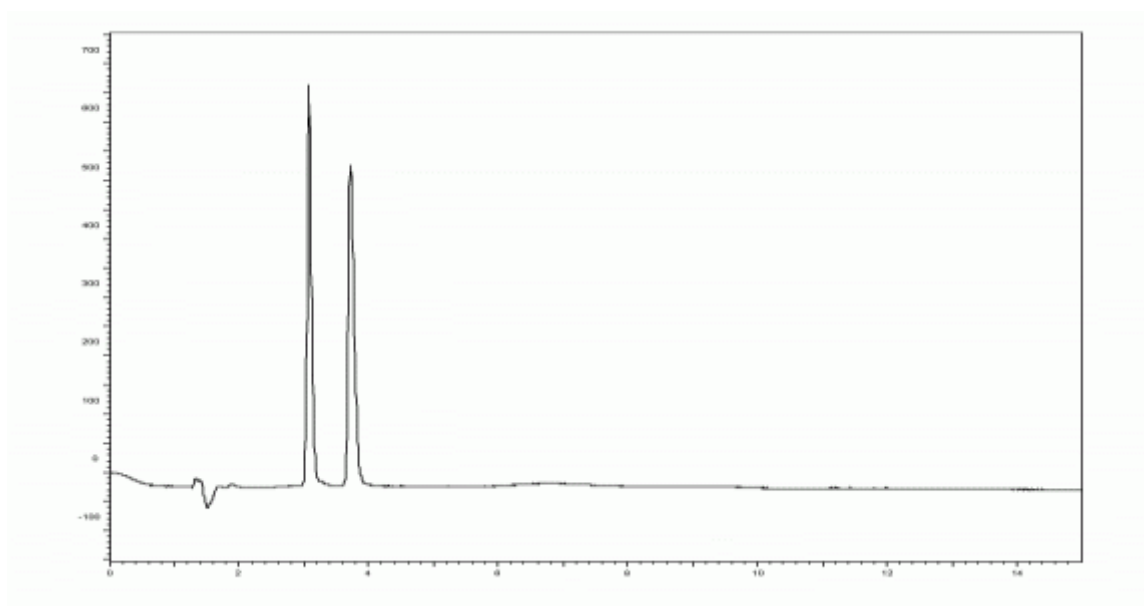

## ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® AY-H 250 x 4.6mm / 5µm
Mobile Phase Composition	CO <sub>2</sub> / 2-propanol = 85 / 15
Chromatographic Mode	SFC

Flow rate (ml/min)	3	Injection amount (mg)	
Temperature (°C)	RT	Pressure (MPa)	12
Detection	UV 210 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	3.07	2.07	1 - 2	4.54
2	3.72	2.72		