



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	3-Methoxy-DL-tyrosine
---------------	-----------------------

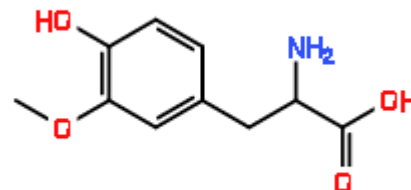

CHROMATOGRAM

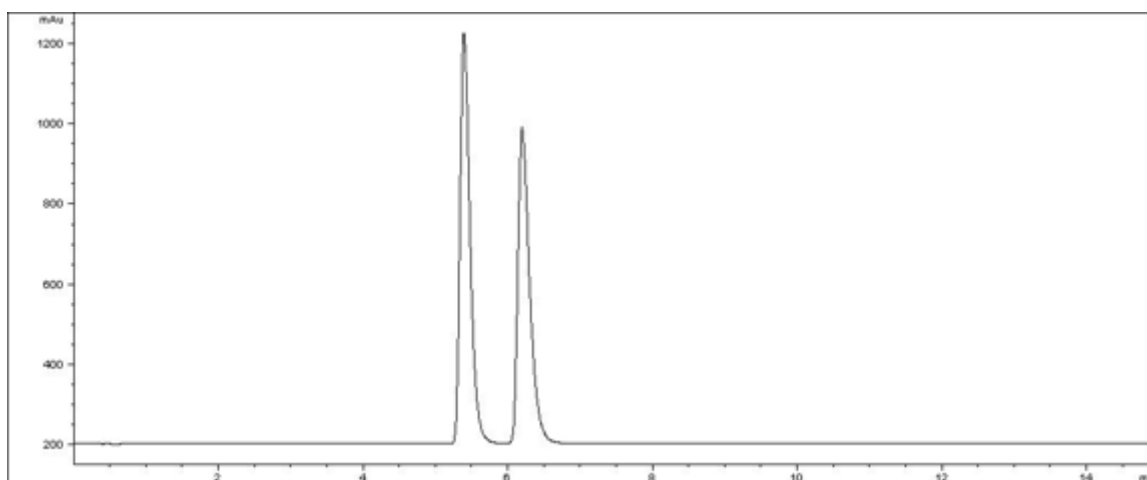

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® ZWIX(+)
Column Size / Particle Size	150 x 3mm / 3µm
Mobile Phase Composition	50mM formic acid + 25mM diethylamine in methanol / acetonitrile / water = 49 / 49 / 2
Chromatographic Mode	Polar Phase

Flow rate (ml/min)	0.5	Injection amount (mg)	0.005
Temperature (°C)	10	Pressure (MPa)	
Detection	ELSD	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	5.39	1.96	1 - 2	3.05
2	6.19	2.4		