



## APPLICATION RESULTS

## COMPOUND IDENTIFICATION

Compound Name	Nipecotnic acid
---------------	-----------------

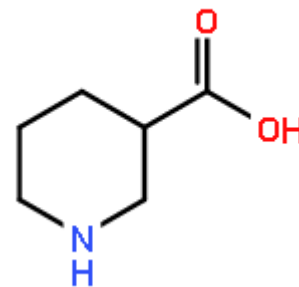

## CHROMATOGRAM

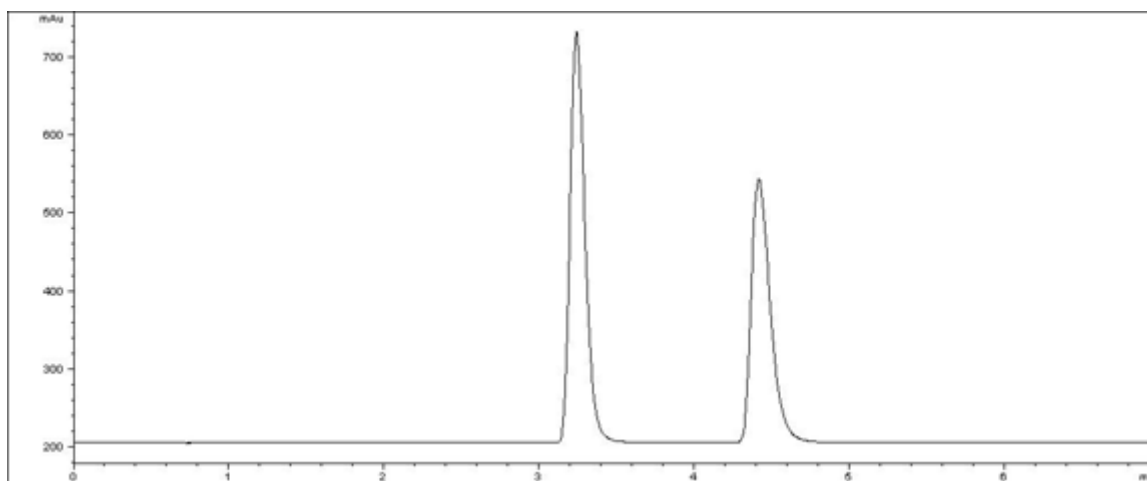

## ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® ZWIX(+) 150 x 3mm / 3µm
Mobile Phase Composition	50mM formic acid + 25mM diethylamine in methanol / water = 90 / 10
Chromatographic Mode	Polar Phase

Flow rate (ml/min)	0.5	Injection amount (mg)	0.0025
Temperature (°C)	25	Pressure (MPa)	
Detection	ELSD	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	3.24	0.79	1 - 2	6.02
2	4.42	1.44		