


APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	Polythiazide
---------------	--------------



CHROMATOGRAM

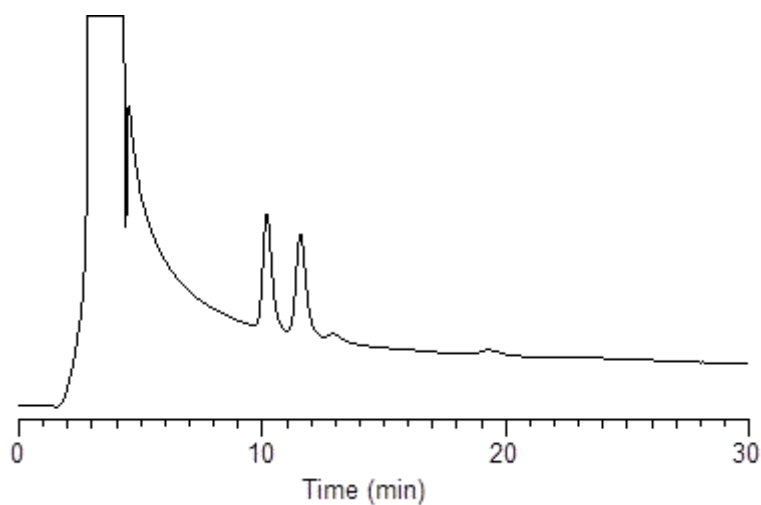

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® IC 250 x 4.6mm / 5µm
Mobile Phase Composition	n-hexane / ethanol / diethylamine = 80 / 20 / 0.1
Chromatographic Mode	Normal Phase

Flow rate (ml/min)	1.0	Injection amount (mg)	0.01
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 230 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	10.2	2.4	1 - 2	1.19
2	11.6	2.9		