



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	Pregabalin
---------------	------------

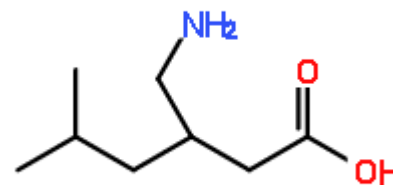

CHROMATOGRAM

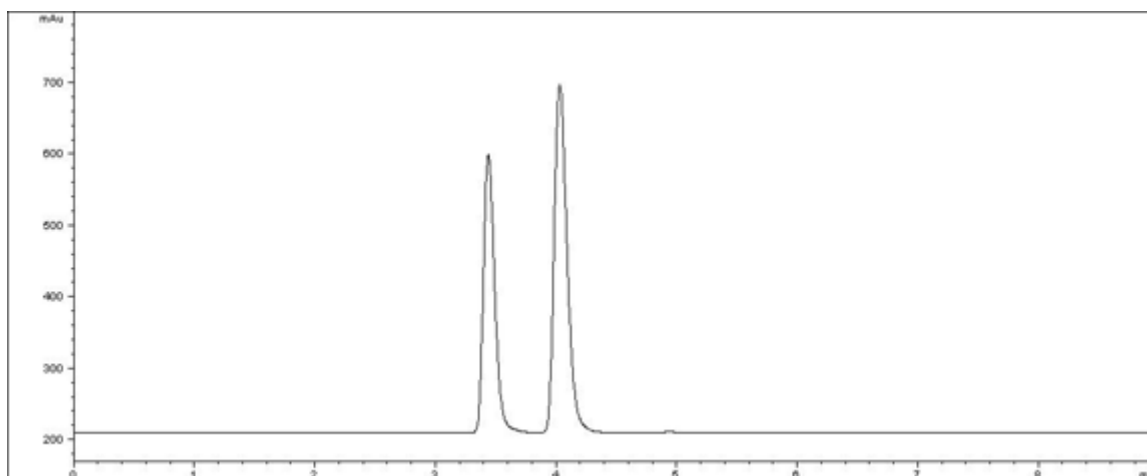

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® ZWIX(+) 150 x 3mm / 3µm
Mobile Phase Composition	50mM formic acid + 25mM diethylamine in methanol / water = 90 / 10
Chromatographic Mode	Polar Phase

Flow rate (ml/min)	0.5	Injection amount (mg)	0.0025
Temperature (°C)	25	Pressure (MPa)	
Detection	ELSD	Elution order	R/S

No.	Rt	K'	Alpha Value	Resolution
1	3.44	0.9	1 - 2	3.32
2	4.03	1.22		