



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	Primaquine
---------------	------------

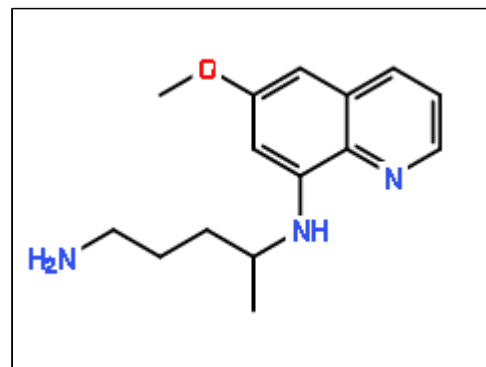

CHROMATOGRAM

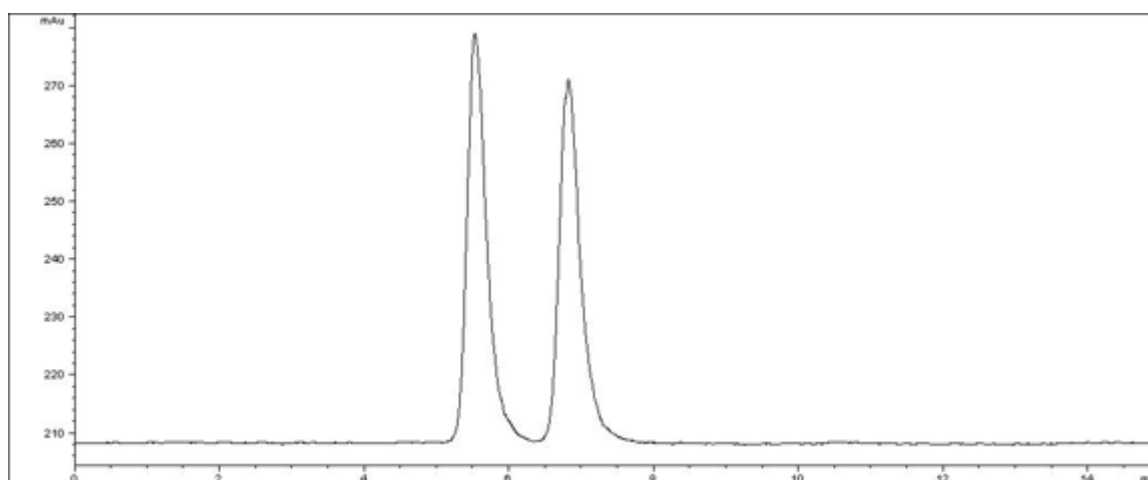

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® ZWIX(+) 150 x 3mm / 3µm
Mobile Phase Composition	12.5mM triethylamine in methanol / water = 98 / 2
Chromatographic Mode	Polar Phase

Flow rate (ml/min)	0.5	Injection amount (mg)	0.0015
Temperature (°C)	25	Pressure (MPa)	
Detection	ELSD	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	5.53	2	1 - 2	2.63
2	6.82	2.71		